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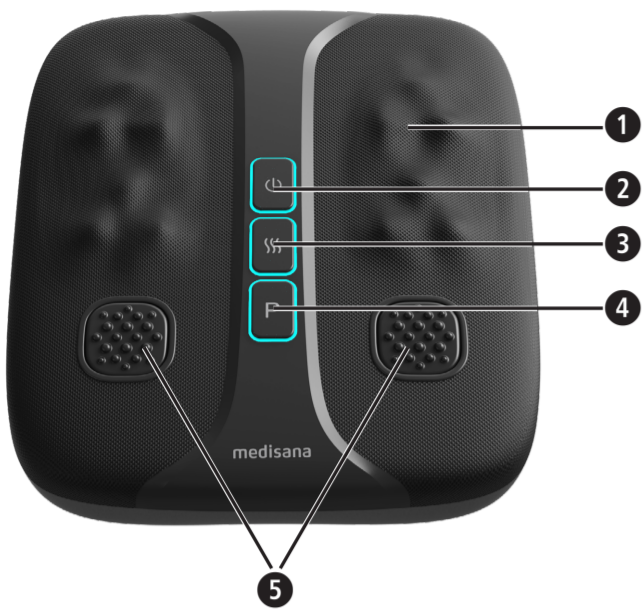
Explanation of symbols

IMPORTANT

Non-observance of these instructions can result in serious injury or damage to the device.

WARNING


These warning notes must be observed to prevent any injury to the user.

CAUTION

These notes must be observed to prevent any damage to the device.

NOTE

These notes give you useful additional information on the installation or operation.

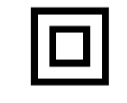Protection category II

LOT

LOT number

Manufacturer

Recycling symbols / codes:

These are used to provide information about the material and its proper use and recycling.

GB Safety Instructions

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

about the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Only insert the mains plug into a socket when the appliance is switched off.
- Keep the cable and the appliance away from heat, hot surfaces, moisture and liquids. Do not touch the mains plug when you are standing in water and always take hold of the plug with dry hands.
- Appliance components carrying voltage may not come into contact with liquid.
- Do not take hold of a device if it falls into water. Pull out the power plug from the socket immediately.
- Connect the unit so that the power plug is freely accessible.
- Always switch off the appliance after use with the ON/OFF-button ② and unplug the power plug.
- To disconnect the unit from the power supply, always pull the plug out of the socket. Do not pull on the cable!
- Do not carry, pull or twist the device on the cable.
- You must not continue to use the unit if the lead or the plug are damaged.
- If the power supply cable is damaged, it may only be replaced by medisana, an authorised dealer or by qualified persons in order to avoid a hazard.
- Ensure that the cable will not cause tripping. The cable must not be kinked, trapped or twisted.

specific groups of people

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision and instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children should be supervised to ensure that they do not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- Do not use the device, or only after consulting your doctor, if:
 - » you are pregnant,
 - » you have a pacemaker, artificial joints or electronic implants,
 - » you suffer from one or more of the following illnesses or complaints: circulatory disorders, varicose veins, open wounds, bruises, skin tears, phlebitis,
- Do not use the appliance near your eyes or near any other sensitive areas.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Take care when using the heating function. The appliance has hot surfaces. Do not use the appliance on small children, vulnerable individuals or people who are not sensitive to heat.
- Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.

before using the unit

- Check the plug, cable and massager carefully for damage before each use. Do not attempt to use a defective unit.
- Do not use the appliance if it, or the power cord show signs of damage, if it is not in perfect working order, or if it has been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.

operating the unit

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using the device.

Only use the massager indoors.

Do not use the massager in damp rooms (e.g. when bathing or showering).

CAUTION

Make sure that the maximum operating time of 15 minutes is not exceeded!

- The device must only be used for its intended purpose as described in the instruction manual.
- Using the unit for any other purpose invalidates the warranty.
- Only use the appliance in the proper way as described in this user manual. The appliance is not designed to support your full body weight! Never stand on the unit!
- Do not use the appliance while you are asleep or when you are lying in bed.
- Do not use the appliance before going to sleep. The massage has a stimulating effect.
- Never use the appliance while driving a vehicle or operating a machine.
- The maximum operating time for one session is 15 minutes. Prolonged use reduces the service life of the appliance.
- Very long periods of use can lead to overheating. So after continuously operating for 15 minutes, user should stop it for at least 15 minutes before using it again.
- Never put the unit down next to an electric stove or other sources of heat.
- Never leave the unit unattended if it is connected to the power supply.
- Do not allow the appliance to come into contact with pointed or sharp objects.
- Do not cover up the appliance when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.

maintenance and cleaning

- The unit is maintenance free. If a fault does occur, just check whether the mains power plug is properly connected.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all guarantee rights.
- Ask your specialist dealer and have repairs carried out by the service only.
- Do not immerse the device in water or other fluids.
- If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.

Scope of delivery

Please check first whether the unit is complete and does not show any damage.

If in doubt, do not operate the unit and contact your dealer or service centre.

The scope of delivery includes:

- 1 **medisana** Foot massager **LM 100**
- 1 Quick Start Guide (short instructions)

If you notice any transport damage when unpacking, please contact your dealer immediately.

Packaging is recyclable or can be returned to the raw material cycle. Please dispose of packaging material that is no longer required properly.

WARNING

Please ensure that the packaging material is kept away from the reach of children! Risk of suffocation!

Use

Place the unit on the floor in front of a comfortable chair so that it cannot tip over and connect the device to your power supply. The front feet of the unit are adjustable if required.

The large buttons on the unit can also be operated with the foot.

- Switch the unit on with the ON/OFF button ②, sit relaxed on the chair and place your feet on the unit on the shiatsu massage heads ①. Place your heels on the vibration plates ⑤. The LED of button ② lights up blue.
- You can switch between the different modes using the mode button ④: Shiatsu massage (standard mode after switching on, LED blue) → Shiatsu and shaking massage slow (after 1st press on mode button ④, LED red) → Shiatsu and shaking massage fast (after 2nd press on mode button ④, LED violet) → Shiatsu massage (standard mode after 3rd press on mode button ④, LED blue).
- Using the heat button ③, the heating function can be switched on in 2 stages: 1st button press → heating function, low level, button lights up red, 2nd button press → heating function, high level, button lights up violet, 3rd button press → heating function off. **The heating function can only be used during a massage mode (not as a single function without massage).**

Switching off the unit

You can switch off the appliance at any time by pressing the ON/OFF button ②. After 15 minutes, the appliance switches off automatically. Then let the appliance cool down for at least 15 minutes before using it again.

Cleaning and care

- Before cleaning the appliance, make sure that the appliance is switched off and the power supply is disconnected. Allow the unit to cool down.
- Clean the massager with a slightly damp sponge only. Never use brushes, aggressive cleaning agents, petrol, thinners or alcohol. Rub the unit dry with a clean dry cloth.
- Never immerse the appliance in water to clean it and ensure water cannot get into the appliance.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.
- Take care to roll up the mains lead to avoid damaging it.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.

Technical specifications

Name and model:	medisana Foot massager LM 100
Power supply:	220-240V~ 50-60Hz 50W
Operating conditions:	Only in dry rooms
Storage conditions:	Clean, dry and cool
Dimensions approx.	368 × 331 × 115 mm
Weight approx.:	3.4 kg
Item number:	88400
EAN number:	40 15588 88400 9

Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH,


Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under <https://docs.medisana.com/88400>.

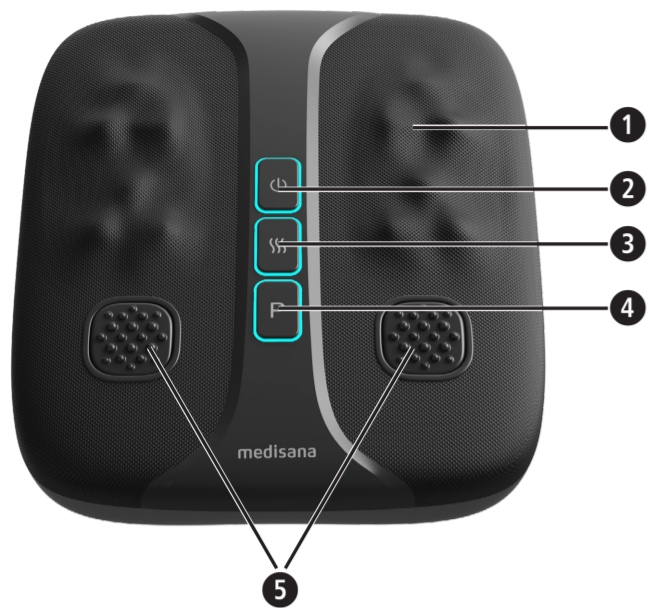Information about service can be found here:
www.medisana.com/servicepartners

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