
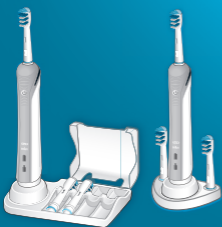
Oral-B®

TRIZONE powered by BRAUN

1000
2000
3000

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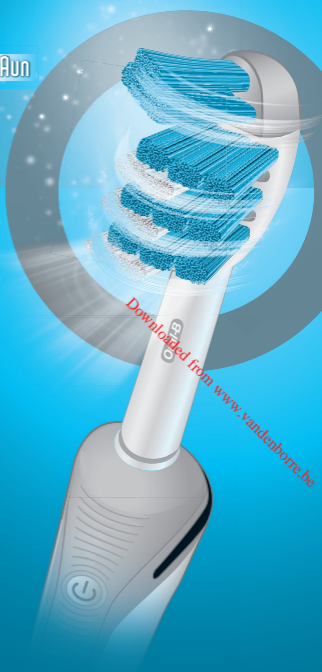English

Thank you for choosing Oral-B TriZone and welcome to Oral-B. This toothbrush has been designed to offer you and your family a unique brushing experience that is both safe and effective.

Before operating this toothbrush, please read these instructions carefully and save this manual for future reference.

IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical or mental capabilities, unless they are supervised by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children.
- Children should be supervised to ensure that they do not play with the appliance.
- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service Centre. Do not disassemble the product except when disposing of the battery.
- When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
- Do not insert any object into any opening of the appliance/charger.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- Do not use on children under age 3. Replace brush (head) every 3 months or sooner if brush head becomes loose due to damage.

Specification

For voltage specifications refer to the bottom of the charging unit.

Description

- a TriZone brush head
- b Pressure sensor light
- c On/off button (mode selection)
- d Handle
- e Charge indicator light
- f Low charge indicator light
- g Basic charger
- h Brush head holder (depending on model)
- i Brush head compartment with protective cover (depending on model)

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the basic charger (g) or the assembled charging unit (see picture B/C) into an electrical outlet and place the toothbrush handle (d) on the charging unit.
- The green charge indicator light (e) flashes while the handle is being charged. Once it is fully charged, the green charge indicator light shines continuously for 5 seconds; then it turns off. A full charge enables up to 7 days of regular brushing (twice a day,

2 minutes) (picture 1).

- When the rechargeable battery is low, the red low charge indicator light (f) starts flashing after turning off the toothbrush.
- For everyday use, the toothbrush handle can be stored on the plugged-in charging unit to maintain it at full power. Overcharging is impossible. However, for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

Please note: If the charge indicator light is not flashing immediately, continue charging. It will appear in 10–15 minutes.

Using your toothbrush

Brushing technique

This brush is designed to work with all common brushing techniques. You can continue to use your current brushing technique, e.g. you can brush just like you would do with a manual toothbrush. You also may consult your dentist or dental hygienist about the right technique for you.

We recommend the following: Wet brush head (a) and apply toothpaste (before switching the appliance on); you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on. You start and stop brushing by pressing the on/off button (c) of the handle (picture 2). Place the toothbrush bristles against the teeth at a slight angle towards the gumline.

Apply light pressure and start brushing in back and forth motions. Brush first the outsides, then the insides. Finally, brush the chewing surfaces with the same back and forth motion throughout the mouth (picture 5). Brush for 2 minutes. Please ensure that you brush all four quadrants of your mouth equally.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

Brushing modes (depending on model)

«Daily Clean»	Exceptional mouth cleaning.
«Sensitive»	Gentle, yet thorough cleaning for sensitive areas.
«Whitening»	Exceptional whitening for occasional or everyday use.

How to switch between modes:

- When pressing the on/off button your toothbrush automatically starts in the «Daily Clean» mode.
- To switch to the «Sensitive» and «Whitening» mode, successively press the on/off button.
- To turn off your toothbrush from the «Whitening» mode, push the on/off button once more.
- If you wish to turn off your toothbrush from the «Daily Clean» or «Sensitive» mode, push and hold down the on/off button until the motor stops.

Professional Timer

A short stuttering sound at 30-second intervals reminds you to brush all four quadrants of your mouth equally (picture 4). A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time. The elapsed brushing time is memorised even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds the timer resets.

Pressure Sensor

Protecting your teeth and gums from hard brushing, your toothbrush has a pressure control feature installed. If too much pressure is applied, the red pressure sensor light (b) will light up reminding you to reduce pressure (picture 3). In addition you will also hear a different sound while brushing. Periodically check the operation of the pressure sensor by pressing lightly on the brush head during

use.

Brush heads

Oral-B TriZone brush head

features triple zone cleaning action for outstanding plaque removal, even in between teeth.

Oral-B also offers you a variety of different Oral-B brush heads that fit your Oral-B toothbrush handle and can be used for tooth-by-tooth precision cleaning.

Oral-B FlossAction brush head

features micropulse bristles that get deep in hard-to-reach areas between your teeth, gently removing hidden plaque.

Oral-B Precision Clean brush head

cups each tooth with a curved bristle trim and reaches deep between teeth at the same time.

Oral-B Sensitive brush head

features an extra soft bristle structure that is gentle on teeth and gums.

Oral-B ProWhite brush head

features a dentist-inspired polishing cup to remove surface stains from coffee, tea and tobacco for effective brightening. Please note that children under 12 years old should not use the Oral-B ProWhite brush head.

Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway down approximately in a 3-month period, signaling it is time to change your brush head. If the bristles splay before the color recedes, you may exert too much pressure on teeth and gums.

We do not recommend using the Oral-B FlossAction or Oral-B ProWhite brush head with braces. You may use the Oral-B Ortho brush head, specifically designed to clean around braces and wires.

Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

Cleaning recommendations

- After brushing, rinse your brush head or accessory carefully under running water with the handle switched on.
- Switch off and remove the brush head or accessory. Clean both

parts separately under running water then wipe them dry before reassembling the toothbrush.

- Disassemble charging unit before cleaning. Never place the charger in dishwasher or water; it should be cleaned with a damp cloth only.
- Brush head holder, brush head compartment and protective cover are dishwasher-safe (picture 6).

Subject to change without notice.

Environmental Notice

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

Guarantee

We grant 1 year guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

Downloaded from www.oral-braun.be