Tefal®

EN

FR

NL

DE

RU

goodvibes[®]
LIFE+

BM962x



Instruction manual - Notice d'utilisation - Gebruiksaanwijzing - Bedienungsanleitung - Manuale di istruzioni - Manual de instrucciones - Manual de instruções - Εγχειρίδιο οδηγιών - Kullanma kılavuzu - Instrukcja obsługi - Uživatelská příručka - Használati útmutató - Návod na použitie - Manual de instrucțiuni - Uputstvo za upotrebu - Uputeza uporabu - Ръководство с инструкции - Руководство по эксплуатации - Priročnik znavidili - Посібник з експлуатації

1

EN - Body composition scale/ FR - Impédancemètre/ NL - Impedantiemeter/
DE - Körperfettanalysewaage/ RU - Весы с измерением жировой мышечной масс

EN - Weighing unit modification can only be done in the app
FR - Le changement d'unités de poids n'est possible qu'à partir de l'application
NL - De weegeenheid kan alleen via de app worden gewijzigd
DE - Die Änderung der Wiegeeinheit kann nur über die App vorgenommen werden
RU - Единицы измерения массы можно настроить только в приложении


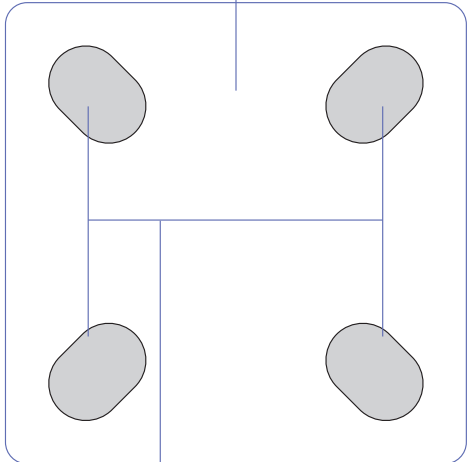
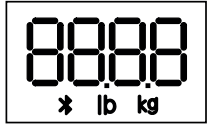
- EN 3 x AAA included
- FR 3 x AAA fournies
- NL Inclusief 3 x AAA
- DE 3 x AAA enthalten
- RU 3 x AAA батарейки входят в комплект

EN Electrodes FR Electrodes NL Electroden DE Elektroden RU Электроды

2

EN - Forbidden/ FR - Interdictions/ NL - Verboden/ DE - Untersagt/ RU - Запрещается

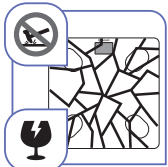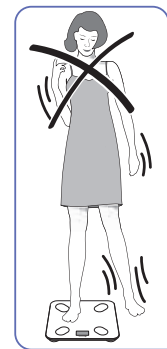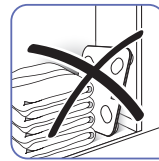EN Thanks for choosing TEFAL Goodvibes body composition Scale. This product employs BIA (bioelectrical impedance analysis) technology for your body composition data and Bluetooth technology to connect the scale to your smartphone. To have accurate and consistent measurement results, please use it at the same time (preferably twice a week, 15 min after waking up) and under same conditions. Please avoid measurements:

- During or immediately after strenuous exercises
- Immediately after a shower, bath or sauna
- With heavy alcohol consumption
- During pregnancy
- During fever

During the measurement, please stand straight and barefoot. Wearing socks or dirty feet can affect the measurement accuracy. The scale should not be used to diagnose or treat any medical condition. All data is for reference and can not be a substitute for advice from a licensed healthcare professional.

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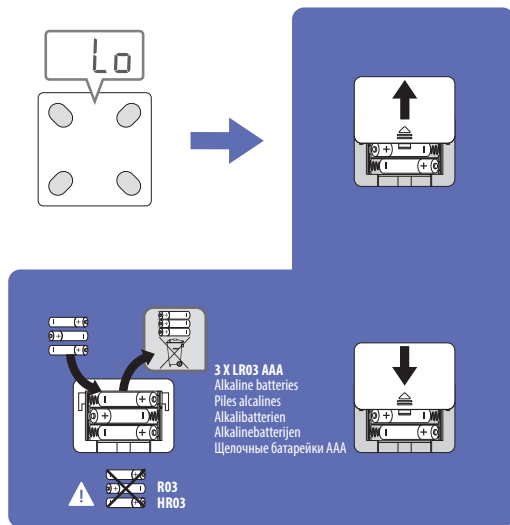
EN WARNING

- Do not attempt to recharge non-rechargeable batteries.
- Do not mix new and used batteries or different types of batteries.
- Ensure that the polarity is correct when inserting the batteries.
- Empty batteries must be removed from the device and disposed of in a safe manner.
- Remove batteries from a device if it is not intended to be used for a long period.
- Do not short-circuit the power supply terminal.

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