GB9212

PAGE

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User's manual facial tanner

Please read these instructions carefully before assembling and using the sun lamp GB9212. Make sure you only use this appliance in accordance with the recommended treatment times, tanning intervals, and tanning distance.

Power supply


This sun lamp should only be connected to a 230V 50Hz socket. Never place the sun lamp in a bathroom, sauna, or any other damp room.

Protective goggles

Ultraviolet radiation is harmful to the eyes. A pair of protective goggles, which must be worn with every use, has therefore been included with this sun lamp. Never use "regular" sunglasses instead of the protective goggles. Even shutting the eyes provides inadequate protection. A spare pair is also included.

Operation

Start the facial tanner by opening the two side panels by pushing button **A** (see drawing **1**). Panels will open themselves by a built in spring. Open the panels till they cannot go any further **1 B**. DON'T FORCE SINCE THAT MAY CAUSE DAMAGE. The panels should stay in a certain angle to give maximum side tanning on the face. Turn the facial in the right position towards the face. The chosen position can be fixed by turning the knob on the right side of the facial (turn to the left **1 E**). Plug the mains lead into the socket and set the timer **1 C** by turning the knob to the right to the time you wish to use the facial tanner.

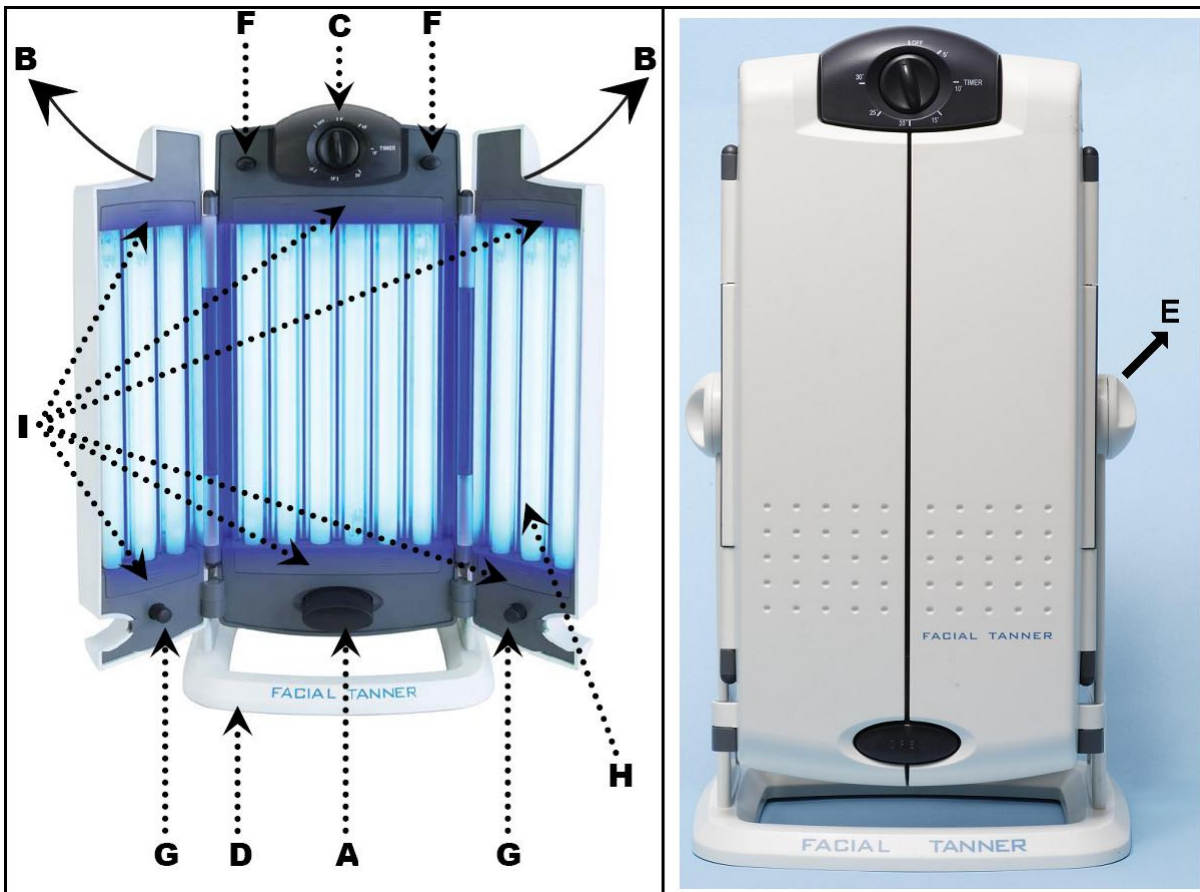DRAWING 1

Timer

Using the timer will automatically switch on the facial and lamps will start burning. Turn the knob to the right to the desired time (max. 30 min.). The timer will then shut off the lamp after the time has elapsed. If you stop using the facial tanner before the time is up just let

the timer run until it ends. DON'T TURN THE TIMER BACK SINCE THAT MAY DAMAGE THE TIMER.

Maintenance, cleaning and technical information

First remove the plug from the socket! Wipe the outside of the sun lamp with a slightly damp cloth. Do not use an abrasive cloth. Make sure that no water gets into the sun lamp. The light tubes can be dusted with a dry cloth or a soft brush. The reflector plate should only be dusted with a very soft brush (if necessary) from time to time. If components such as the lamps or power cord (if damaged) need to be replaced the work should always be carried out by a qualified electrician to avoid dangerous situations. The lamps have a very long operating life of about 300 hours. This facial tanner has been built according to all European safety standards, which allows it to carry the CE-mark. Insulation class II, UV A radiation type 3, mains voltage 230 Volts, 50 Hz, power 200 Watts.

General information about tanning

You can tan safely all year round underneath this sun lamp if you follow these guidelines correctly. Our skin is tanned by invisible ultraviolet radiation from the sun. This radiation darkens the pigment already present in our skin and stimulates the formation of more pigment known as melanin. The production of new melanin cells and the darkening of these make the skin stay tanned longer. There are a number of different types of ultraviolet radiation. Ultraviolet A radiation (UVA), which colors the pigment cells brown, and ultraviolet B radiation (UVB), which stimulates the production of melanin cells, but also causes sunburn, are the most well-known types. Your sun lamp provides mostly UVA radiation so that the danger of sunburn is precluded. How quickly and how deep you tan depends on your skin type. The skin type also detain lines how long one can (safely) stay in the sun. This is true for both the real sun as well as the sun lamp. In order to make the best use of your sun lamp, first check your skin type using the following table. Further on in these instructions, you will find the different skin types in the tables with the recommended tanning times. If your skin falls into the type 1 category, we recommend you do not use the facial tanner at all.

SKINTYPE

Skin type Characteristics	Type 1	Type 2	Type 3	Type 4
Eye color?	Light blue, grey or green	Blue, grey or green	light brown	Dark brown
Hair color?	Red, blond	blond	Dark blonde, brown	Brown, black
Skin color?	Reddish, pale	pale	Pale, light brown	Brown
Freckles	lots	some	A little	None
Tanning experience	Doesn't or barely tans	Tans lightly	Tans reasonably	Tans easily
Sunburns	Burns quickly and peels	Peels quickly	Occasionally burns and doesn't peel	Rarely or never burns

The tanning program for facial Eurosolar GB9212

Before you begin, read all the instructions!

SKINTYPE 2		SKINTYPE 3		SKINTYPE 4	
Treatment	Treatment in minutes	treatment	Treatment in minutes	treatment	Treatment in minutes
1	10	1	10	1	10
2	10	2	10	2	10
3	10	3	10	3	15
4	10	4	15	4	15
5	15	5	15	5	15
6	15	6	20	6	20
7	15	7	20	7	25
8	15	8	20	8	30
9	20	9	25	9	30
10	20	10	25	10	30
11	20	11	30	11	30
12	20	12	30	12	30

Between two treatments there should be at least 48 hours!

Maximum number of treatments:

- Skin type 2 60
- Skin type 3 43
- Skin type 4 33

These recommended hours only apply for use with the recommended light tubes. Determine your skin type using the table. The maximum number of minutes for one tanning session can be found there.

Keep a minimum of 15 cm distance between your face and the lamp. A greater distance is allowed although the tanning effect will diminish drastically. Always use the protective goggles! Should your skin become red and/or tight and sensitive, it is recommended you decrease the treatment time by 25%. After a treatment of 12 days, a maintenance tanning time of 30 minutes a week is sufficient to keep your tan but don't exceed the maximum number of treatments. Unplug the lead after each use. If you wish you can close the panels by pushing them slightly. The panels will be locked automatically. BEFORE CLOSING THE PANELS LAMPS SHOULD BE TURNED OFF!. If this by accident is forgotten the two safety switches will automatically shut down the operation. However don't use this as a regular procedure!

WARNINGS!!

- Never use the facial tanning lamp in the bathroom, sauna, or any other damp room.
- If any of the following health conditions apply to you, we suggest you consult a doctor before using the facial tanning lamp: diabetes, epilepsy, pregnancy, prescription of medicine or drugs, UV hypersensitivity, prone to dizziness, prone to fainting, high blood pressure, skin disease.
- If your skin forms large freckles, if birthmarks get larger, or if in doubt about certain changes in your skin, we advise that you consult a doctor. People (having been) afflicted with skin cancer or those with a disposition towards contracting it should not use the facial tanner.
- Certain other conditions medication, drugs, and cosmetics can make you more sensitive to UV radiation, if in doubt, consult your doctor before use of the sun lamp. Do not use this medication when using the facial tanner and thoroughly remove all cosmetics before use. Do not use tanning creams or oils when using the facial tanner.

- Do not use the facial-tanning lamp if your skin is still red from sunburn.
- Always take a one-day break after each tanning session. Do not sunbathe outdoors on the same day you have used the facial tanner.
- UV radiation can damage the skin and eyes. Therefore, be sure to follow these instructions closely and always use the protective goggles provided.
- Use the facial tanning lamp alone in a room. People without protective goggles should stay away from the facial tanning lamp while it is being used.
- Ultraviolet radiation from the sun or UV lamps can damage skin and eyes. These biological effects depend on the quality and quantity of the radiation and on the individual sensitivity of the eyes and skin.
- Excessive exposure to UV radiation can cause sunburn, age the skin prematurely, and increase the risk of skin tumors.
- Never use the facial tanning lamp if the timer, mains lead, or other electrical mechanisms are damaged or do not function properly. Consult the dealer where you bought the product.
- Replace the light tubes with the recommended type Cosmedico Cosmolux X 15W/10043. Instructions how to change the lamps please read chapter "changing lamps" further on. The reflector panel should never be replaced by another type.
- Keep objects subject to discoloration out of the facial tanning lamp's radiation, such as curtains, bed linen, furniture, paintings, etc. If necessary, cover these with a sheet.
- Never let the facial tanner be used by children or those who, due to a handicap, cannot operate the device themselves and/or those who cannot follow the recommendations in these instructions.

CHANGING THE LAMPS


If tanning results are decreasing a lot (mostly after 300 hours use) you may change the lamps. ONLY USE **COSMOLUX X 15W/10043**

1. make sure the facial tanner is switched off and the power cord has been taken out of the socket
2. remove the covers (**I** on drawing **1**) by pushing them in the direction indicated on the plastic
3. remove the tubes by turning them slightly and carefully
4. push in the new tubes and fix them by turning them into opposite direction
5. put the covers back in place

Parts and functions see drawing 1

- A Release button to open the panels
- C Timer with knob
- D Stand

- E Turning point with fixation knobs
- F Safety switches (turning off the lamps when the panels are accidentally closed during operation. (NO STANDARD PROCEDURE). See instruction manual.
- G Spring buttons to release the panels easier
- H Lamps 15W Cosmolux X
- I Lamp covers

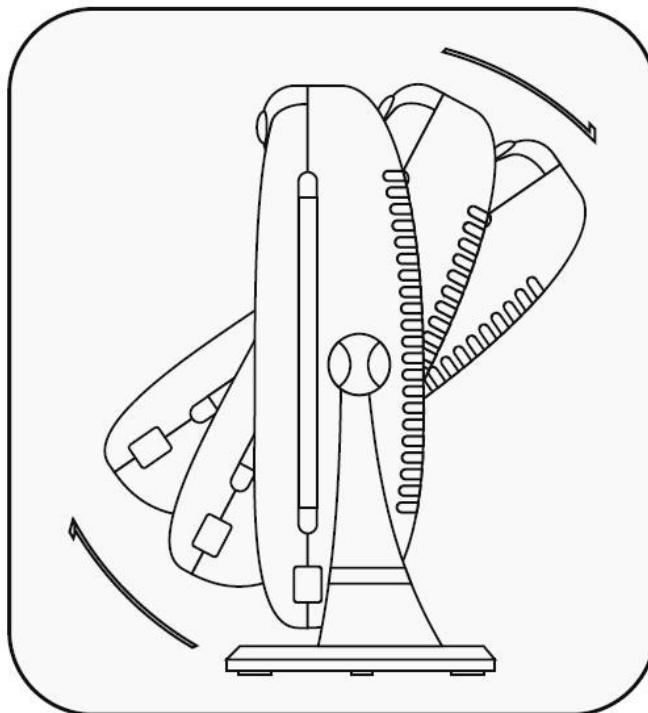Different positions o adjust with fixation knob E

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PRESAT S.A.