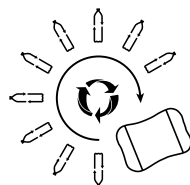
LANAFORM

COMFORT

Hottle

ELECTRIC RECHARGEABLE HOT-WATER BOTTLE

Outer layer
made out
of recycled
PET material

EN Charging time: $\pm 10'$
FR Temps de chargement: ± 10
NL Laadtijd: $\pm 10'$
DE Ladedauer: $\pm 10'$

EN Up to 3hrs use
FR Jusqu'à 3h d'utilisation
NL Gebruiksduur tot 3u
DE Bis zu 3h Einsatz

EN Protection against overheating
FR Protection contre la surchauffe
NL Bescherming tegen oververhitting
DE Überhitzungsschutz

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1

ENGLISH

Thank you for choosing the LANAFORM electric hand warmer. Convenient and easy to use, it will keep you warm for several hours by charging it in just a few minutes. Furthermore, you can use it completely safely thanks to its superior protection.

Finally, the cover of this hot-water bottle is partially made out of fabric, the fibres of which are obtained from recycled PET material, from plastic bottles.

▲ Please read all the instructions before using your electric hot water bottle, especially this basic safety information.

☑ The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.

Always follow these instructions when using this appliance.

Check that your mains voltage is the same as that of your device.

This appliance is not intended for use by people (including children) with limited physical, sensory or mental abilities, or by people lacking experience or knowledge, unless they are supervised or have been instructed beforehand on how to use the appliance by someone responsible for their safety. Children should be supervised to make sure they do not play with the device.

Do not use accessories other than those recommended by LANAFORM or accessories which are not supplied with this device.

The electrical hot water bottle must be used with the charger supplied in its packaging.

Do not connect or disconnect the plug of the device with wet hands.

If the power cable is damaged, it must be replaced by a special lead or similar lead available from the supplier or its after-sales department.

Check the electrical hot water bottle regularly. Do not use this device if the plug is damaged, if it is not working properly, if it has been dropped on the floor, or if it is damaged or has fallen into water. In this event, have the device examined and repaired by the supplier or its after-sales department.

If it falls into water, unplug it immediately before retrieving it.

Check that the electrical hot water bottle is fully sealed before using it.

Do not use the electrical hot water bottle when it is plugged in.

The hot water bottle must be charged on a flat surface and not be vertical or sloping.

Do not use this device under a cover or under a cushion. Excessive heat can cause fire, electrocution or injury.

Take all necessary precautions to avoid burns.

Handle the electrical hot water bottle with care. The internal temperature may reach in fact 75°C and 60°C on its surface once charged. If necessary, you may wait until the temperature drops before using the electrical hot water bottle.

Please remove the charger and do not use the electrical hot water bottle if you hear an unusual noise or notice an unusual smell as well as a leak from the electrical hot water bottle.

To make sure the heat is evenly spread, you are advised to shake the electrical hot water bottle gently after 5-6 minutes when charging it.

Electrical appliances should never be left unattended when connected. Unplug it when you are not using it.

Keep the power cable away from hot surfaces.

Always unplug the device after use or before cleaning it.

Do not insert any object into the electrical hot water bottle. This could impair its operation.

Do not hit or cut the surface of the hot water bottle. Do not exert too much pressure on the electrical hot water bottle when you are using it as this may cause the device to leak a leak or damage to it.

Disconnect the electrical hot water bottle if it overheats.

If the electrical hot water bottle does not automatically stop charging after 15 minutes (the light remains on), disconnect the charger plug.

Do not move the electrical hot-water bottle using the charger's power cable.

Before storing the electrical hot water bottle, make sure its temperature has dropped back down.

Do not crush the electrical hot water bottle by putting items on top of it once it is been stored.

The electrical hot water bottle is intended for domestic use only. Do not use it in a hospital environment or for any kind of medical treatment.

The electrical hot water bottle is not intended to be used by children.

Do not use the electrical hot water bottle on people who are sick, on babies or on those insensitive to heat.

Do not use the electrical hot water bottle to warm up your feet. You must not walk on the electrical hot water bottle, crush it or exert pressure on it.

Do not use the electrical hot water bottle in a microwave or conventional oven.

You are advised not to let the electrical hot water bottle come into direct contact with your skin.

If, when using the electrical hot water bottle, you notice unusual red marks appearing on your skin, or even discomfort or pain, stop using it for at least an hour. If the problem persists, do not continue using the hot water bottle and, if necessary, consult your doctor immediately.

Never use this appliance in a room in which aerosols (sprays) are used or in a room where oxygen is administered.

FEATURES

- The body of the electrical hot water bottle is pliable and very comfortable.
- The walls consist of 2 thin, thermo-compressed film layers, increasing its durability and safety.
- With the technology it has, the electrical hot water bottle has a faster heating system, therefore heat is maintained for longer.

INSTRUCTIONS FOR USE

When you have removed the product from its packaging, check the integrity of the device and that of the components supplied with it. If in doubt, do not use the product and contact the product's dealer. Do not use the product if there are signs of defects and return it to your dealer.

Charging the electrical hot water bottle

- 1 Lay the electrical hot water bottle down flat on a table.
- 2 Insert it into the charger as shown on the diagram (■ 1). Connect the charger plug to the socket. The light comes on and the electrical hot water bottle heats up.
- 3 Once the electrical hot water bottle reaches an internal temperature of 75°C, the light goes out and you can disconnect the charger. You can then use the electrical hot water bottle.

❶ If you want the electrical hot water bottle to stay hot for longer, squeeze the electrical hot water bottle to mix the water, wait for a few minutes and then continue charging it.

❷ It takes about 10 minutes to charge the electrical hot water bottle (depending on the ambient temperature) so that the surface of the electrical hot water bottle reaches a temperature of $\pm 60^{\circ}\text{C}$.

MAINTENANCE AND CLEANING

Storage

Put the electrical hot water bottle in its box or in a dry place.

Avoid contact with sharp edges and other pointed objects, which might cut or pierce the fabric surface.

To prevent the power cable breaking at any point, do not wind it round the latter.

Make sure that the electrical hot water bottle is completely dry before storing it

Cleaning

Unplug the electrical hot water bottle from the charger plug and wait a few minutes before cleaning it.

Clean it using a soft, slightly damp sponge.

Do not immerse it in liquid to clean it.

Never use abrasive detergents, brushes or a thinning agent to clean it.

TECHNICAL FEATURES

- Voltage: 220-240V
- Frequency: 50 Hz
- Power: 440W
- Composition of the outer lining: 25% recycled polyester, 75% cotton
- Composition of the inner lining: 100% polypropylene

⚠️ BREAKDOWNS AND FIXES

Ⓢ The hot-water bottle inflates and does not heat up properly.

- Ⓢ Changes in atmospheric pressure may cause air to accumulate in the hot-water bottle.
- Ⓢ Carefully remove excess air from the hot-water bottle by removing the plug located behind the connector cover (unscrew).

🗑️ ADVICE REGARDING THE DISPOSAL OF WASTE

All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.

When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

LIMITED WARRANTY

LANAFORM guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM's control.

LANAFORM may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM's Service Centre. Any attempt to maintain this product by a person other than LANAFORM's Service Centre will render this guarantee void.

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