LANAFORM®

Luminescence

Light care

Downloaded from www.vandenborre.be

Light effects with
health benefits

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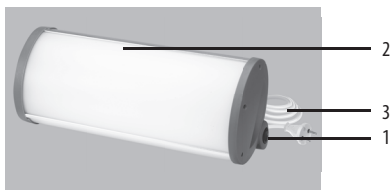
fig. 1

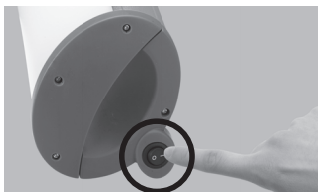
fig. 2

fig. 3

fig. 4

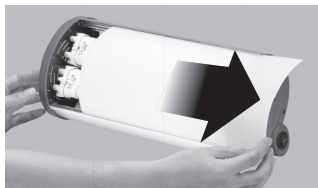
fig. 5

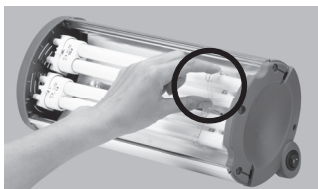
fig. 6

fig. 7

Thank you for purchasing the LANAFORM® LUMINESCENCE lamp. The LUMINESCENCE lamp works by using a very clear light source that simulates natural sunlight.

Thanks to the technology and the significant luminous intensity (approx. 10,000 lux at 17 cm)*, the LUMINESCENCE lamp will help quickly restore your energy and vitality.

You can use the LUMINESCENCE while carrying out other activities such as reading, writing, using the telephone, or even working.

* Refer to the usage instructions for more detailed information.

ABOUT LIGHT EXPOSURE

Light and its intensity have a major influence on our bodies. The sun's rays are weaker in autumn and winter. Traditional artificial lighting in buildings is incapable of compensating for this deficit in natural light.

Numerous healthcare professionals agree on the benefits to the body of exposure to light. By light, they are referring to daylight of a certain intensity. In winter, light intensity does not exceed 2500 lux*, whereas it can reach 100,000 lux on a clear day in summer.

LANAFORM has developed this light-emitting device specifically to enable you to benefit from an additional dosage of light during the autumn and winter months.

Exposure to light consists of using a light with a strong lux coefficient* to simulate natural sunlight.

* The lux is the unit of measurement for the luminance of light

IMPORTANT SAFETY INFORMATION

Please read these instructions before using the device. Please also observe these instructions carefully when using the device.

- Only use this device for the type of use described in this manual.
- This device is not intended for use by persons (including children) with reduced physical, sensory or mental faculties, or persons lacking in experience or knowledge, unless they are supervised or given prior instructions on using the device by a person responsible for their safety. Particular attention must be taken around children to ensure that they do not play with the device.
- Before plugging in the device, check that the voltage indicated on the device corresponds to the voltage of the local network.
- Do not use accessories that are not recommended by LANAFORM® or that are not supplied with this device.
- No modifications may be made to this device.
- If the power cable has been damaged, it must be replaced with a similar one available from the supplier or the supplier's aftersales service.
- Do not use this device if the plug has been damaged, if it is not working properly, if it has fallen onto the floor, if it has been damaged or if it has fallen into water. Have the device checked and

- repaired by the supplier or the supplier's aftersales service.
- Do not carry this device by the power cable and do not use the cable as a handle.
 - Always unplug the device after use or before cleaning it.
 - Electrical equipment must never remain plugged in when unsupervised. Unplug it when you are not using it.
 - Keep the electric cable away from sources of heat.
 - Do not use this product before going to bed. Exposure to light can have a stimulating effect that prevents sleep.
 - Never let anything fall into any of the openings, and never insert an object into them.
 - Do not use this device in a damp environment (in a bathroom, near a shower etc).
 - Never use this device in a room where aerosol products (sprays) are used or in a room where oxygen is administered.
 - Do not use this device under a cover or cushion. An excess of heat may cause a fire, electrocution or injury to the user.
 - This device is designed purely for use in the home. Use of the LUMINESCENCE is not equivalent to medical treatment.
 - Do not use the device in a room that is not illuminated by another source of light - this may lead to ocular fatigue. Use this device in tandem with normal lighting (daylight or artificial light).
 - If you feel any pain while using this device, stop using it immediately and consult your doctor.
 - We advise you not to look at the lamp during your session. You may experience sore eyes during the initial sessions; this will soon disappear.
 - Consult your doctor before starting to use the LUMINESCENCE lamp if:
 - you suffer from or have suffered from chronic depression
 - you suffer from hypersensitivity to light
 - you have an eye problem
 - you are taking certain medications (e.g. certain antidepressants, sedatives, malaria tablets) that might have an adverse effect on use of the LUMINESCENCE lamp
 - you have high blood pressure.
 - The period of exposure essentially depends on the distance between the LUMINESCENCE lamp and you. Refer to the section entitled «Usage instructions» to work out the recommended distance and period of exposure.

COMPONENTS OF THE PRODUCT

- 1) On/Off button
 0. Off
 1. 2 x 24W low energy fluorescent PLL tubes.
- 2) Reflective screen
- 3) Electric cable

USAGE INSTRUCTIONS

- For best results, it is recommended that the LUMINESCENCE lamp be used between the hours of 6 a.m. and 8 p.m.

Installation:

1. Place the LUMINESCENCE on a solid surface such as a table.
2. Plug the power cable into the power socket.

Usage :

- Press the On/Off button (Fig. 2) to illuminate the LUMINESCENCE lamp.
- The recommended period of use for the LUMINESCENCE depends on the distance to the lamp. Nevertheless, it is generally recommended that the device be used for 2 hours at a luminous intensity of 2500 lux. You can use the table below to easily calculate the recommended period of exposure:

Distance	Position 1 (2 tubes)	Recommended period
+/- 17 cm	10.000	30 mins
+/- 25 cm	5.000	1 hr
+/- 45 cm	2.500	2 hrs

- Position yourself in front of the device.
 - o The distance between your face and the screen must be between 17 and 50 cm. Adjust this to the distance that feels best for you and for your eyes in particular.
 - o You are not required to sit directly opposite the device.
 - o You can split the duration of your luminothérapie session.
- Do not look directly at the lamp. During the session, you may read, write, use the telephone etc.
- The exposure must be regular over a prolonged period. During the darkest days of the year, repeat use for at least 7 consecutive days, if not more, depending on your personal requirements.
- After use, stop the lamp by switching the «On/Off» button located on the top of the device to

position «0» and unplug it.

- If you feel eye or other pain during the light exposure session, stop using the LUMINESCENCE lamp immediately and consult your doctor.

MAINTENANCE AND STORAGE

Cleaning:

- Unplug the device and let it cool down.
- Use a cloth or damp sponge to clean the exterior surface of the lamp.
- Never use a scrubbing brush, abrasive products or aggressive detergents such as acetone or white spirit to clean the device.

REPLACING THE TUBES

The device contains no parts needing replacement other than the tubes. The tubes are expected to last for approximately 8000 hours of use (which corresponds to several years of normal use). The tubes must be replaced by similar ones.

To replace the tubes, please follow the following instructions:

- Using a screwdriver (not supplied), release the 3 screws located at the back of the product (see Fig. 3).
- Then remove the reflector screen as illustrated in Figs. 4 & 5.
- Then, pull back the fixing clip and remove it as shown in Figs. 6 & 7.
- Once you have replaced the tube, replace the reflective screen before using the LUMINESCENCE lamp.

TECHNICAL SPECIFICATIONS

Voltage: 220-240V

Frequency: 50Hz

Power consumption: 48W

TROUBLESHOOTING

Problem	Possible cause	Solution
No light	No current	The device may not be plugged in correctly or your power source may be defective
	Main switch is off	Switch on main switch
One side of the product is not working	Defective neon tube	Change the neon tube
	Neon tube at end of life	Change the neon tube
The light continues to flicker for 20 seconds after switching on	Incorrectly positioned neon tube	Reposition the neon tube

ADVICE REGARDING THE DISPOSAL OF WASTE

- All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.
- When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

LIMITED WARRANTY

LANAFORM® guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM® guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM® product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM®'s control.

LANAFORM® may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM® will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM®'s Service Centre®. Any attempt to maintain this product by a person other than LANAFORM®'s Service Centre will render this guarantee void.

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