LANAFORM®

Baby Stone

Children's cosy corner

Version 1 du 01/03/2016

Figure 1

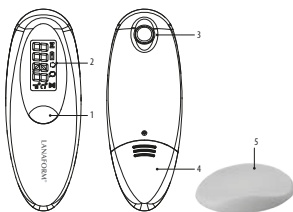
Figure 2

Figure 3

Figure 4

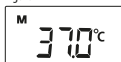
Figure 5

Figure 6

Figure 7


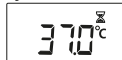
Figure 8

EN : TITLE AND DESCRIPTION OF DISPLAY SYMBOLS ON the brochure
describing the product's features:

Symbol Definition	Identification	Information provided by the manufacturer
Manufacturer's details		Geon Corporation _ No. 12, Gung Ye Road, Hsi Hu Chang Hwa Hsien - Taiwan
Batch number		
BF device		
Instructions for use		Follow the instructions on how to use the device
Protective rating liquid and solid particle	IP 22	A coding system to indicate the degrees of protection provided by an enclosure against access to hazardous parts, ingress of solid foreign objects, ingress of water and to give additional information in connection with such protection.
Please note: see instructions for use		Please note: consult the enclosed documents. Please read these instructions before using the device. Please also follow these instructions carefully when using the device.
Notified body		
European representative		MYM Soporte tecnico y de Calidad, S.L. Avda. de los Rosales, 32, 28935, Mostoles (Madrid), Spain
Electrical and electronic waste (WEEE)		The device should not be disposed of with other household waste at the end of its serviceable life.

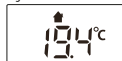PLEASE READ ALL THE INSTRUCTIONS BEFORE USING YOUR THERMOMETER, ESPECIALLY THIS BASIC SAFETY INFORMATION:

- When using this product, please make sure you follow all the instructions listed below. Any deviation from these instructions may cause injury or affect the accuracy of the readings.
- Check that your mains voltage is the same as that of your device.
- Always keep the monitor away from heat sources such as stoves, radiators, etc.
- Do not use this device in a room in which aerosols (sprays) are used or in a room

in which oxygen is administered.

- Do not place the thermometer in direct sunlight, or expose it to high temperatures, humidity or dirt.
- Do not immerse the device in water.
- Do not drop the device and do not use it if it is damaged.
- This device must be used at an ambient temperature of between 15°C et 35°C.
- Do not expose the thermometer to electric shocks.
- Do not expose this thermometer to extremes of temperature above 55°C or below -25°C.
- Do not use the device if the relative humidity is above 95% or below 10%.
- Do not use the device in the vicinity of large electro-magnetic fields, such as cordless or portable telephones.
- The accuracy of the measuring process may be disrupted if the forehead is covered by hair, a cap or a scarf.
- Please wait at least 30 minutes after taking a shower, exercising, eating or drinking before taking any temperature.
- It is recommended that the thermometer should be kept at room temperature for around 10 minutes before use.
- Use of this thermometer is not a substitute for consulting your doctor.
- If there is a problem with your device, please contact your reseller.
- Do not try to repair the device yourself.
- Do not take the device as a whole apart.
- This device is not intended for use by people, including children, with limited physical, sensory or mental abilities, or by people lacking experience or knowledge, unless they are supervised or have been instructed beforehand on how to use the equipment by someone responsible for their safety. Children should be supervised to make sure they do not play with the device.
- Only use the thermometer for its intended purpose, as described in these instructions.
- The thermometer is pre-set when it leaves the factory. You do not therefore have to calibrate the device when you first use it.
- Put the thermometer back on its base and keep out of the reach of children after it has been used.
- The forehead thermometer has a probe which is used to take the temperature on the forehead only. Do not use it to take a temperature in the rectum, mouth, ear or armpit.
- Body temperature is generally lower in the morning and gradually rises throughout the day. Please therefore take your body temperature at the same time during the day and under the same conditions of use.
- The temperature must be taken between the eyebrows and roots of the hair.
- Do not place the thermometer on a scar or any problem areas of the skin, open wounds or abrasions or if you are suffering from any physical trauma.
- Consult your doctor if you are taking any medicines. These might increase your body temperature.
- The device's performance may be adversely affected if the thermometer is used or stored outdoors.
- The thermometer should only be used indoors.

FEATURES

SNAPSHOT : Short waiting time: temperature taken in degrees Celsius just under 5 seconds.

COMFORTABLE : Its ergonomic shape ensures it is comfortable to use.

ECONOMICAL : The thermometer's standby mode has a battery life of at least 1000 usages.

DURABLE : Probe detection technology provides the system with a long serviceable life.

ACCURATE : The temperature is accurate to $\pm 0.1^{\circ}\text{C}$.

COMPACT : The thermometer's compact, lightweight design makes it the ideal tool, both at home and when travelling.

USER-FRIENDLY : The data that appear on the liquid crystal screen are easy to understand. Backlit screen for visibility at night.

INSTRUCTIONS FOR USE (Figure 1)

1 : ON/OFF button 2 : Backlit LCD screen 3 : Probe 4 : Battery compartment
5 : Base

HOW TO USE THE THERMOMETER

MEANING OF THE THERMOMETER'S SYMBOLS AND SIGNALS

Temperature settings:  Celsius: $^{\circ}\text{C}$ - Fahrenheit: $^{\circ}\text{F}$. Memorising the last recorded result: «  » . Reading the temperature taken on the forehead: «  »

Thermometer's ready-to-read temperature: «  » . Continuous display for 10 seconds of the temperature taken: «  »

Low-battery alarm:  . Standby mode and continuous display of ambient temperature: «  »

OPERATING THE THERMOMETER

The "BABY STONE" forehead thermometer by LANAFORM® is a high-tech electronic device for taking body temperature. Please remove the protective film covering the battery compartment before using your thermometer for the first time. Press the thermometer's on/off button to switch it on (Figure 2) All the existing icons will appear on the screen for a second (Figure 3). The last temperature taken is then displayed on the screen showing «  » (memory) for a second (Figure 4).

A beep tone followed by the «  » icon tells you that you can now take your child's temperature (Figure 5) .To do this, carefully place the thermometer's probe on your baby's forehead as shown in Figure 6 and keep it there for 5 seconds.

An audible alert will sound and the temperature as well as an hourglass symbol «  » will be instantly displayed on the screen for ten seconds or so (Figure 7). If the child's temperature rises above 37.8°C (100°F), 6 consecutive audible alerts will sound. You then have the option of taking the baby's temperature again immediately since a beep tone followed by the icon «  » will again be displayed on the screen.

After one minute's inactivity, the thermometer will automatically switch to standby mode, displaying the ambient temperature of the room and the icon of a house «  » will be displayed on the screen (Figure 8). If you want to take your child's temperature again, press the on/off button to activate the mode for taking the forehead temperature. To switch off the thermometer, press this on/off button again. To choose the $^{\circ}\text{C}$ or $^{\circ}\text{F}$ unit of measurement, switch the thermometer off. Then keep your finger on the thermometer's on/off button for 5 seconds. The $^{\circ}\text{C}$ or $^{\circ}\text{F}$ unit of measurement will be displayed on the screen. Select the temperature required by pressing this same on/off button. Let the thermometer to set itself once you have chosen the unit of measurement and this

will be displayed on the screen. Do not press the on/off button as you might then change the unit of measurement you have chosen.

REPLACING THE BATTERY

1. Press the cover protecting the battery and slide it towards the outer edge before removing it. 2. Get a new lithium CR2032 battery (3V) ready. 3. Remove the old battery from where it is located and replace it with the new one. 4. Replace the protective cover.

 -To avoid accidents and misuse, keep the batteries out of the reach of children. - Please dispose of used batteries in the appropriate waste collection bins. To prevent environmental pollution, do not simply discard batteries. - Do not recharge the batteries. - Do not swallow the batteries. - Do not expose the batteries to fire or they might explode. - Do not dismantle the device except for its battery compartment. - It is recommended that the batteries be removed when the thermometer is not being used for long periods.

The following two tables are a guide showing average human temperatures according to the place where the temperatures are taken and age.

Range of normal temperatures according to the measuring method used.

Measuring method	Normal temperatures in °C
In the rectum	36.6 to 38
In the mouth	35.5 to 37.5
Under the armpit	34.7 to 37.3
In the ear	35.8 to 38
On the forehead	35.8 to 37.8

Normal temperatures according to age

Age	Temp. °C
0-2 years	36.4 to 38.0
3 to 10 years	36.1 to 37.8
11 to 65 years	35.9 to 37.6
>65 years	35.8 to 37.5

Notes: 1. The temperature of a body in good health varies according to a lot of factors: age, time of day, physical activity, outside temperature, etc., and more specifically when taken on the forehead (influenced by type of skin (greasy or dry, wrinkled or otherwise), ambient temperature (air conditioning, etc.). It is therefore important to establish the range of usual and normal temperatures for each member of your family to provide benchmarks for detecting any fever.

2. During any consultation, it is important to let your doctor know that the temperature taken was on the forehead.

TECHNICAL FEATURES

Normal conditions of use	
Display resolution	0.1 °C
Operating temperature	15 to 35°C (59°F to 95°F) / 10% - 95% RH / 700 - 1060hpa
Storage temperature	-25 to 55°C (-13°F to 131°F) / 30% - 85% RH / 700 - 1060hpa

(Battery) power supply	3V continuous (1 «CR2032» lithium battery)
Battery life	+/- 1000 measurements (in normal use)
Level of protection against solid particles and liquids	IP 22
The thermometer is not water-proof	

Measuring range	35 to 43°C
Measuring time	5 seconds
Area measured	Area where temperature is taken is between the eyebrows and roots of the hair.
Accuracy of measurement	$\pm 0.2^{\circ}\text{C}$ ($< 35.5^{\circ}\text{C}$) $\pm 0.1^{\circ}\text{C}$ ($35.5^{\circ}\text{C} \rightarrow 42^{\circ}\text{C}$) $\pm 0.2^{\circ}\text{C}$ ($42.1^{\circ}\text{C} \rightarrow 43^{\circ}\text{C}$)
Memory	Last measurement displayed when switched on

MAINTENANCE AND CLEANING

- When you start taking a new measurement and in order to ensure accuracy and good hygiene, always clean the tip of the probe with a cotton bud moistened with alcohol. Wait for two minutes before taking the temperature again.
- Use a soft, dry cloth to clean the body of the thermometer. Never use a rough towel and never immerse the thermometer in water or other liquids.

TROUBLESHOOTING

Anomalies	Possible reasons	Solutions
1. The screen constantly displays the temperature gauge as too low "Lo"	a) Make sure that the tip of the probe is clean b) The temperature taken is below 35°C (95°F)	a) Clean the probe with a soft, dry cloth or cotton bud moistened with alcohol. b) Wait until the "OK" symbol appears and start the operation again
2. The screen constantly displays the temperature gauge as too high "Hi"	a) Make sure that the tip of the probe is clean b) The temperature taken is more than 43°C (109°F)	a) Clean the probe with a soft, dry cloth or cotton bud moistened with alcohol b) Wait until the "OK" symbol appears and start the operation again
3. The low-battery symbol «  » appears on the screen	The battery is low	Replace it with a new battery
4. The screen constantly shows "Err 1"	The temperature shown is beyond the thermometer's detection limit	-Start the operation again once the "OK" icon appears on the screen -Put the thermometer in the room for 30 minutes and then start the operation again

5.The screen constantly shows "Err 3"	There is a problem with the thermometer and it is no longer able to function as it should	Contact your reseller
---------------------------------------	---	-----------------------

Advice regarding the disposal of waste:

- All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.
- When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law. Prior to disposal, please remove the battery and dispose of it at a collection point so that it can be recycled. Please note: under no circumstances should spent batteries be disposed of with household waste!!!

LIMITED WARRANTY

LANAFORM® guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below. The LANAFORM® guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM® product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM®'s control. LANAFORM® may not be held liable for any type of circumstantial, indirect or specific damage. All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied. On receipt, LANAFORM® will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM®'s Service Centre®. Any attempt to maintain this product by a person other than LANAFORM®'s Service Centre will render this guarantee void.

The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.

DECLARATION DE CONFORMITY CE 0197
Medical device - Category IIa.
We, the Geon Corporation, hereby declare that the medical device «Baby Stone» LA090111 complies with European Directive 93/42/EC concerning medical devices and its most recent amendment 2007/47/EC.

CONTACT INFORMATION

The device is manufactured by: Geon Corporation N°12, Gung Ye Road, Hsi Hu Chang Hwa Hsien, Taiwan

The device is distributed by: Lanaform SA Rue de La légende 55, B-4141 Belgium
 Authorised European representative: MYM Soporte tecnico y de Calidad, S.L.
 Avda. de los Rosales, 32, 28935, Mostoles (Madrid), Spain

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Geon Corporation N°12, Gung Ye Road, Hsi Hu Chang Hwa Hsien, Taiwan Prístroj distribuuje spoločnosť: