PR-S90 connect

Gebrauchsanweisung
Instruction Manual
Gebruiksaanwijzing
Mode d'emploi
Instrucciones de manejo
Οδηγίες χρήσης
Istruzioni per l'uso

**Bitte sorgfältig lesen!
Please read carefully!
A.u.b. zorgvuldig lezen!
A lire attentivement s.v.p.!**
**¡Por favor lea con cuidado!
Παρακαλούμε διαβάσετε με προσοχή!
Da leggere con attenzione!**

CE

MEDISANA®

Art. 88810

PR-S90 connect

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Device and controls

Toestel en bedieningselementen

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Tel.: +41 091 980 49 72 • Fax: +41 091 605 37 55

E-Mail: info@bluepoint-service.ch • Internet: www.bluepoint-service.ch

IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

1 Safety Instructions !

- These scales are not suitable for people with pacemakers or other medical implants. The body fat information may be inaccurate in people with diabetes or other medical/physical limitations. The same also applies to very highly trained athletes.
- This device is not suitable for pregnant women!
- The scale has been manufactured for domestic use. It is not suitable for commercial use in hospitals or other medical institutions.
- Do not use the scale in a location where extremely high temperatures or humidity may occur. Keep the scale away from water.
- Do not place the scale on an uneven surface. Do not stand on the corners/edges of the scale. When weighing yourself, do not stand on one side or the edge of the scales.
- Do not stand on the scale with wet feet. Do not stand on the scale wearing socks.
- The scales have a measuring range up to 180 kg, 396 lb, 28 st. Do not overload the scales – this can cause irreversible damage.
- Handle the scale with care. The scale must not be subjected to impacts or vibration. Do not fall or drop anything on the scales.
- Do not attempt to repair the device yourself in the event of a malfunction since this will invalidate the warranty. Repairs should only be carried out by authorised service centres.
- Children shall not play with the appliance.
- **WARNING!** Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!
- Keep **batteries** out of children's reach! Do not disassemble **batteries**! Do not attempt to recharge these **batteries**! Do not short circuit! Do not throw into a fire.
If a battery has been swallowed seek medical attention immediately!

2 Items supplied and packaging

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use it and contact your supplier or your service centre.

The following parts are included:

- 1 **MEDISANA** body analysis scale **PR-S90 connect**
- 4 x 1.5 V batteries (type AAA)
- 1 Instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required. If you notice any transport damage during unpacking, please contact your supplier without delay.

WARNING Please ensure that the polythene packing is kept away from the reach of children!
Risk of suffocation!

3 Insert/replace batteries

Before you start up your digital scale, insert the four provided 1,5 V batteries (type AAA) into the device. To do this, open the battery compartment **1** beneath the device and insert the batteries. Make sure the polarity is correct. Replace the batterie if the battery change symbol "**LO**" appears in the display or if nothing appears in the display **2** after you have stepped on the device.

Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

Please remove the battery before disposing of the device/unit. Do not dispose of old batteries with your household waste, but at a battery collection station at a recycling site or in a shop.

4 Change the unit of weight

You may have your weight displayed in kilograms (kg), pounds (lb) or stones (st). Turn on the scale. When display show "**0.0**", press **▲ 3**- or **▼ 5** button to change unit KG-LB-ST:LB.

5 Weighing "Step-on"-function

1. Make sure the scales are on a firm, level surface. If you only want to weigh yourself but have moved the scales just before, you must first initialise the scales. To do this, press your foot on the middle of the surface of the scales. "**0.0 kg**" appears in the display **2**. **When the scales switch off, they are ready for the "step-on" function.** This procedure is not necessary if you have not moved the scales.
2. Step onto the scales and stand still. The device switches on automatically and the display **2** briefly shows "**0.0 kg**".
3. Your weight is measured, flashes three times and is then finally displayed.
4. Step down from the scales. The scales switch off automatically after approx.10 seconds.

6 Programming personal data

The scales can weigh, measure and save (except for P0 mode) data (gender, size, age) for 8 people + guest mode (P0). The scales have preset values.

These are: **kg, female, 25 years old, 165 cm.**

1. Press the SET button **4** to switch on the scales. The memory location is flashing in the display
2. Now press **▲ 3** - or **▼ 5** - button to select the desired memory location (**0 - 8**) for the user profile.
3. Press the SET-button to confirm the selection. Now the symbol for the sex is flashing in the display.
 = male  = female  = male/athlete  = female/athlete
4. Set your sex by pressing the **▲ 3** - or **▼ 5** - button to select the appropriate symbol in the display.
5. Press the SET button to save the setting. The preset height starts flashing in the display.
6. Select your height using the **▲ 3** - or **▼ 5** - buttons.
7. Press the SET button to save the setting. The preset age starts flashing in the display.
8. Select your age using the using the **▲ 3** - or **▼ 5** - buttons.
9. Press the SET button to save the setting.
10. The settings are now concluded. In the display "0.0" is shown.

7 Weighing and measuring

1. Press the SET button **4** to switch on the scales. The memory location flashes in the display.
2. By pressing the DOWN- and UP - buttons you may choose the memory location (**0 - 8**).
3. Your personal data is displayed.
4. Wait until "0.0" appears on the display. Step barefoot onto the scales and stand still. Place your feet on the electrodes **6**. First, your weight is displayed.
5. The display shows the values measured for the body fat and body water in one display, and then KCAL (BMR), muscle mass, bone mass value will be shown in one display after. The series of measured values are repeated 3 times in the display. The scales then switch off automatically, step off the scales.

NOTE:

Auto user selection by your last weight. Scale auto chooses user once standing barefoot on the platform, then measuring the body values.

Body fat proportion range in %

| Age | Female (%) | | | | Male (%) | | | |
|--------------|------------|-----------|-----------|-----------|----------|-----------|-----------|-----------|
| | low | normal | high | very high | gering | normal | hoch | sehr hoch |
| 6-15 | 3.0-16.0 | 16.1-21.0 | 21.1-29.6 | 29.7-50 | 3.0-13.0 | 13.1-17.4 | 17.5-25.6 | 25.7-50 |
| 16-30 | 3.0-16.4 | 16.5-22.0 | 22.1-30.4 | 30.5-50 | 3.0-13.6 | 13.7-18.0 | 18.1-26.4 | 26.5-50 |
| 31-60 | 3.0-17.0 | 17.1-22.4 | 22.5-31.4 | 31.5-50 | 3.0-14.0 | 14.1-18.6 | 18.7-27.0 | 27.1-50 |
| 61-80 | 3.0-17.6 | 17.7-23.0 | 23.1-32 | 32.1-50 | 3.0-14.4 | 14.5-19.0 | 19.1-28 | 28.1-50 |

NOTE:

The data in the table is for orientation purposes only. If you have any questions about your proportion of body fat, please consult your doctor.

8 VitaDock+ App*

Weighing yourself and transferring data via Bluetooth to the VitaDock+ app

Make the personal settings on the scales as described in „6 Programming personal settings“. To transfer the data, the scales and the smartphone must be connected via Bluetooth® 4.0 or Bluetooth® Smart®. To do this, install the latest version of the VitaDock+ app.

1. Free download of the VitaDock+ app in the app store or Google Play store.

***System requirements:**

iOS: iPhone 4S and later, iPad 3 and later.




Android: Version 4.3 (Jelly Bean)

Bluetooth® 4.0 technology.

2. Activate Bluetooth® on your iOS or Android appliance.
3. Start the app. Ensure that there is a stable connection to the internet (3G or WLAN). Register now with your VitaDock online account. If you don't have your own account yet, you can register free now at VitaDock Online.
4. Touch the Bluetooth symbol to open the list of devices.

5. Scroll down the device list until you see the scales.

6. Touch the arrow beside the scales to start the setup process.

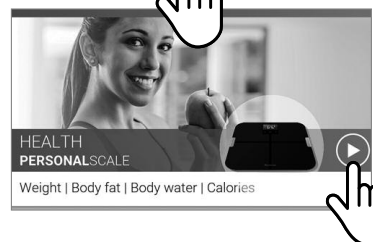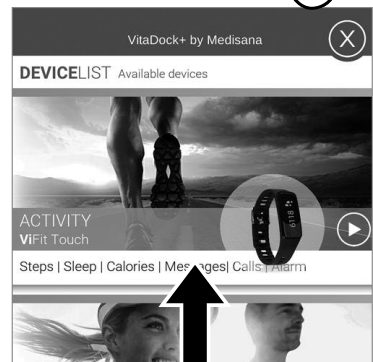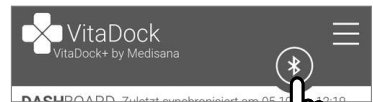7. Press the “Set up scales” button to start the setup process.







8. Then scroll through to choose your target weight. Tap “OK” to confirm the set target weight.

9. You are then prompted to set up a user profile. Tap the ? symbol if you need more information. Touch “OK” to start the connection process

10. You are prompted to weigh yourself. Stand barefoot on the scales. After you have weighed yourself, the scales automatically activate the Bluetooth connection. The app starts synchronisation automatically as soon as it is connected to the scales.

11. The synchronisation symbol rotates while the data is being transferred.

12. Once synchronisation is complete, you can add a note to the latest weighing. Scroll down and select “OK” to save the additional note.

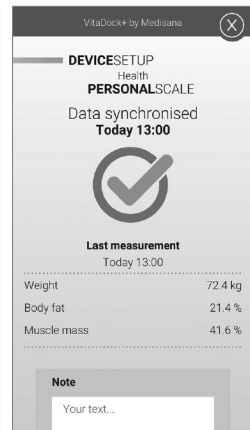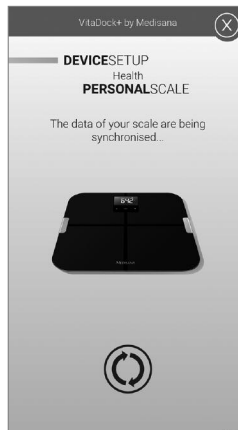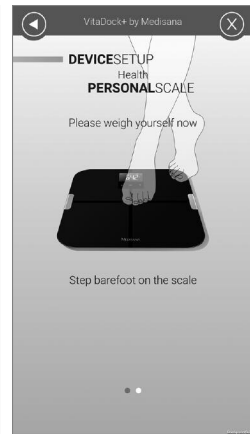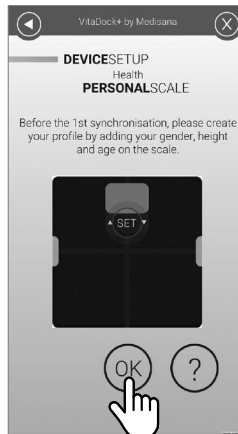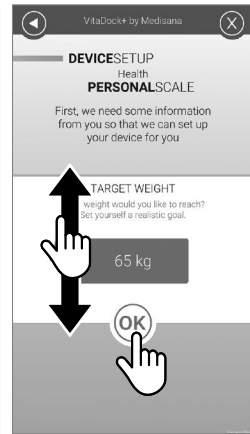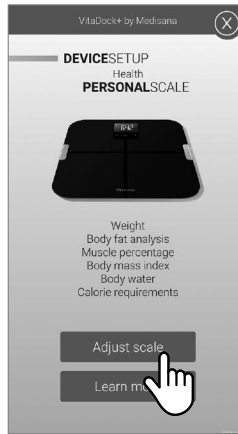In the dashboard, you can touch the graphs or the arrow for the detailed overview. This shows how your weight is progressing. The star marks your set target weight.

Scroll down to see the list of the most recent weighings. You can delete individual weighings or add notes to them.

NOTE: If you use the scales in guest mode „PO“, no data is saved. The weight data cannot be transferred by Bluetooth in this mode.

In the course of continuous product optimisation, we reserve the right to regularly update the VitaDock+ app. Therefore, we recommend that you periodically perform an update. The included images serve only as a reference.

9 Error messages

0-LD = The scale is overloaded.

LO = The battery is low and must be replaced.

If your personal scales do not function to your satisfaction, check the following before contacting a service centre:

- Check whether the batteries have been inserted correctly.
- Check whether the correct unit of weight has been selected.
- Check whether the scales are freestanding on a solid and even surface. They must not be touching a wall or any other object.
- Weigh yourself again.

10 Care and maintenance

Never use aggressive cleaning agents or stiff brushes. Clean the scale using a soft and lightly moistened cloth. Never use abrasive cleaning agents or alcohol. Do not allow any water to get into the unit. Do not use the unit again until it is completely dry.

11 Technical specifications

| | |
|-----------------------------|--|
| Name and model | : MEDISANA body analysis scale PR-S90 connect |
| Voltage supply | : 6 V [≠] , 4 x 1.5 V batteries (type AAA) |
| Display system | : Digital display |
| Memory | : for 8 persons + guest mode (P0). |
| Measuring range | : up to 180 kg, 396 lb or 28 st |
| Maximum measuring deviation | : ± 1 % |
| Body fat | : 3 - 50 % in steps of 0.1 % |
| Body water | : 20 - 75 % in steps of 0.1 % |
| Muscle mass | : 20 - 75 % in steps of 0.1 % |
| Bone mass value | : 0.5 - 8 kg in steps of 0.1 kg |
| Graduation | : 100 g, 0.2 lb or 1/4 st |
| Autom. switch-off | : after approx. 10 seconds |
| Dimensions (L x W x H) | : 319 x 319 x 23 mm |
| Weight approx. | : 2,0 kg |
| Item number | : 88810 |
| EAN number | : 40 15588 88810 6 |

Guidelines / Standards

We hereby declare that bathroom scales model PR-S90 art. 88810 complies with the fundamental requirements of the European RTTE Directive 1999/5/EC. You can request the complete EC Conformity Declaration from Medisana AG, Jagenbergstrasse 19, 41468 Neuss, Germany, or you can also download it from the Medisana homepage.

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under www.medisana.com

12 Warranty and repair terms

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for MEDISANA products is three years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be corrected free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. nonobservance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear (batteries etc.).
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

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MEDISANA AG