BRAUN

Oral-B

D10

Type 4 733

Kids' Power Toothbrush






Braun Infoline

Deutsch	3, 16, 19
English	4, 16, 19
Français	5, 16, 19
Español	6, 16, 19
Português	7, 16, 19
Italiano	8, 17, 19
Nederlands	9, 17, 19
Dansk	10, 17, 19
Norsk	11, 17, 19
Svenska	12, 17, 19
Suomi	14, 17, 19


Internet:

www.braun.com
www.oralb.com


- D** Haben Sie Fragen zu diesem Produkt?
- A** Rufen Sie an (gebührenfrei):
00 800 27 28 64 63
- CH** **08 44 - 88 40 10**
- E** Servicio al consumidor para España:
901 11 61 84
- P** Serviço ao Consumidor para Portugal:
808 20 00 33
- I** Servizio consumatori:
02-6 67 86 23
- NL** Heeft u vragen over dit produkt?
Bel Braun/Oral-B Consumenten-infolijn:
0800-4 45 53 88
- B** Vous avez des questions sur ce produit?
Appelez Braun Belgique
02-711 92 11
- DK** Har du spørgsmål om produktet? Så ring
70 15 00 13
- N** Spørsmål om dette produktet? Ring
22-72-88-10
- S** Frågor om apparaten? Ring Kundservice
020 - 21 33 21
- FIN** Onko Sinulla kysyttävää tuotteesta? Soita
09 452 872

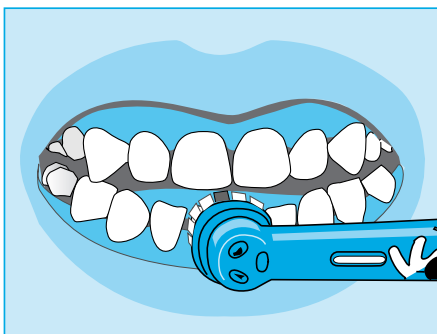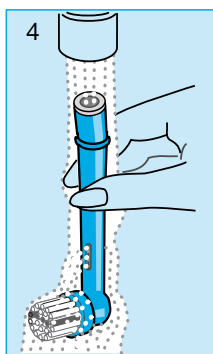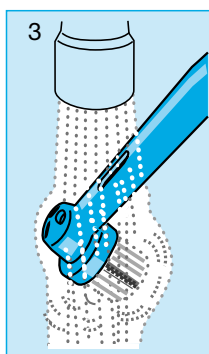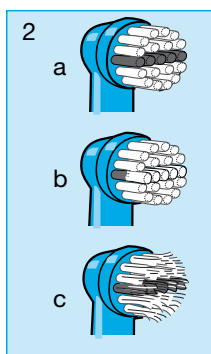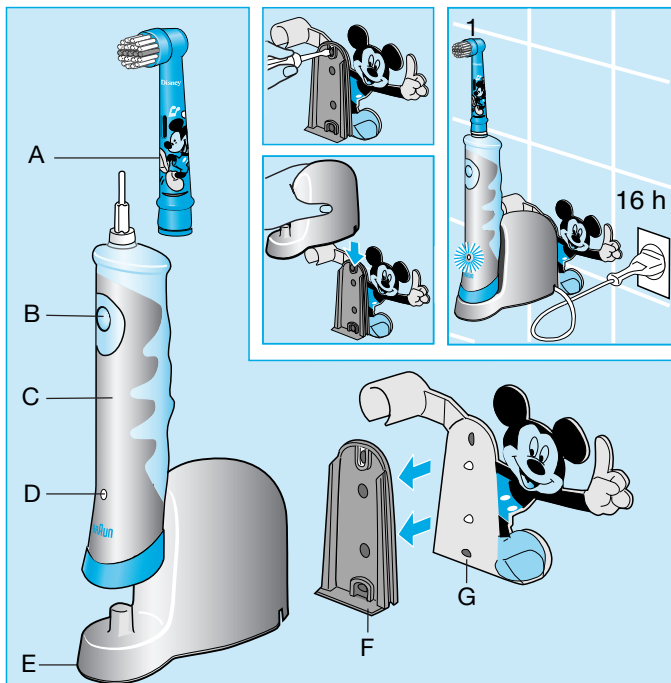outsides

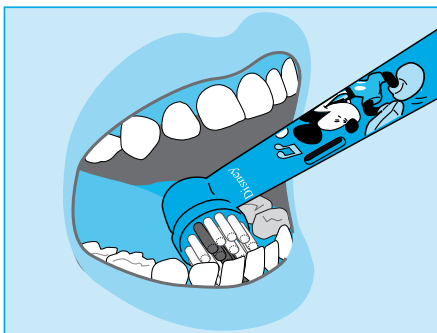insides

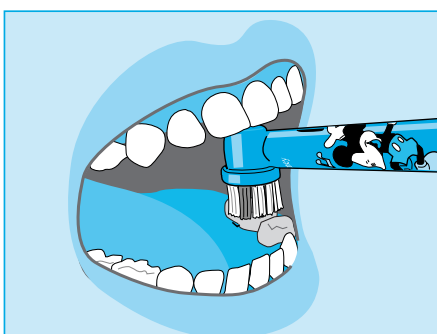tops

Deutsch

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English

Dear parents,
The Braun Oral-B Kids' PowerToothbrush has been developed with leading dental professionals to effectively remove plaque from children's teeth. It will make it easier for your child to brush properly and thoroughly. We hope your child will enjoy his/her new Braun Oral-B appliance.

This appliance is suitable for children 3 years or older. Please read the use instructions carefully and completely before showing your child how to use the appliance. We recommend that parents supervise their children when they start using the appliance in order to help them become familiar with it.

Important

Periodically check the cord for damage. Should this be the case, take the charging unit (E) to a Braun Service Centre. A damaged or non-functioning unit should no longer be used.

Description

- A Brushhead
- B On/off switch
- C Squish Grip handle
- D Charge indicator
- E Charging unit
- F Wall mount
- G Brushhead holder

Specifications

Power supply: see bottom of the charging unit
Handle voltage: 1.2 V

Connecting and charging

The plaque remover is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit. The charge indicator will light up (1).
- A full charge takes 16 hours and will provide an operation time of approximately one hour.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the appliance at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

Removing plaque with a music timer

Use the following steps to show and explain to your child how to use the appliance properly:

- When the handle is taken from the plugged-in charging unit, a short melody will sound.
- Put a small amount of toothpaste onto the brushhead. (We recommend fluoride toothpaste.)
- First guide the brushhead to the teeth, only then turn on the handle with switch (B).
- Start with the lower teeth: Slowly guide the brushhead from tooth to tooth;

- first the outsides,
- then the insides,
- finally the chewing surfaces.

Do not scrub or press too hard. The Braun Oral-B Kids' Power Toothbrush automatically provides the proper cleaning action.

- After one minute of operation, the brush movement will be interrupted by a short melody. This tells your child to move on to the upper teeth. Clean them – again tooth by tooth – until the next melody plays (after another minute of operation). This way you can make sure that your child brushes his/her teeth for at least 2 minutes as recommended by dental professionals.
- Turn off the handle with switch (B) while the brushhead is still in the mouth.

We recommend that you replace the brushhead every three months. The brushhead is provided with a row of blue INDICATOR[®] bristles (2a). With proper brushing twice a day for two minutes, the colour will disappear half-way down these bristles in a 3-month period (2b), provided toothpaste is used. If bristles begin to splay before the colour starts to recede, this indicates that your child exerts too much pressure on teeth and gums (2c).

Replacement brushheads EB10 Available at your retailer or at Braun Service Centres.

Cleaning

After use, rinse the brushhead for several seconds under running water with the handle switched on (3). Then switch off the handle and take off the brushhead. Rinse both parts separately under running water (4) and wipe them dry before putting the handle back onto the charging unit.

Environmental notice

This product is provided with a nickel-hydrate rechargeable battery which does not contain any harmful heavy metals. Nevertheless, in the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Centre.

However, if you prefer to dispose of the battery yourself, please refer to the chapter «Battery removal» on page 15.

Subject to change without notice.

This product conforms to the European Directives EMC 89/336/EEC and Low Voltage 73/23/EEC.

At the end of the product's useful life, please dispose of it at appropriate collection points provided in your country.

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English

Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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Great Britain

Gillette Group UK Ltd.,
Braun Consumer Service,
Aylesbury Road,
Thame OX9 3AX
Oxfordshire
☎ 1800 783 70 10

Greece

Berson S.A.,
47, Agamemnonos,
17675 Kallithea Athens,
☎ (210)-9 47 87 00

Guadeloupe

Ets. André Haan S.A.,
Zone Industrielle
B.P. 335,
97161 Pointe-à-Pitre,
☎ 26 68 48

Hong Kong

Audio Supplies Company,
Room 506,
St. George Building,
2 ICE House Street,
Hong Kong,
☎ (852) 2524 9377

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Tehran Bouran Company
145 Ghaem Magham Farahani
Tehran
☎ +982 18 31 27 66

Ireland (Republic of)

Gavin's Electronics,
83-84, Lower Camden Street,
Dublin 2,
☎ 1800 509 448

Israel

S. Schestowitz Ltd.,
8 Shacham Str.,
Tel-Aviv, 49517,
☎ 1 800 335 959

Italia

Servizio Consumatori Braun
Gillette Group Italy S.p.A.,
Via G.B. Pirelli, 18,
20124 Milano,
☎ 02 / 6678623

Jordan

Interbrands
Wadi Al-Sir, Adwiat alhekma
St. Mekadabi bld.
Amman
☎ +692 582 75 67

Kenya

Radbone-Clark Kenya Ltd.,
P.O. Box 40833,
Mombasa Road,
Nairobi,
☎ 2 82 12 76

Korea

Gillette Korea Ltd.
144-27 Samsung-dong,
Kangnam-ku,
Seoul, Korea,
☎ 080-920-6000

Kuwait

Union Trading Company,
Braun Service Center,
P.O. Box 28 Safat,
Safat Code 13001, Kuwait,
☎ 04 83 32 74

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Philippines

Gillette Philippines Inc.,
Corporate Corner Commerce Avenues
20/F Tower 1, IL Corporate Centre
1770 Muntinlupa city
☎ 02-771071 02-06/-16

Poland

Gillette Group Poland
Sp. z o.o.
ul. Domaniewska 41,
02-672 Warszawa,
☎ 22 548 89 74

Portugal

Grupo Gillette Portugal, Lda.,
Braun Service,
Rua Tomás da Fonseca,
Torre G-2ºA,
1600-209 Lisboa,
☎ 808 2 000 33

Réunion

Dindar Confort,
Rte du Gymnase,
Boite Postale 278,
97940 St. Clotilde,
☎ 026 92 32 03

Romania

Gillette Romania srl.
Calea Floreasca nr. 133-137
et 1, sect 1,
714011 Bucuresti
☎ 01-2319656

Russia

RTC Sovinservice,
Rusakovskaya 7,
107140 Moscow,
☎ (095) 264 41 61

Saudi Arabia

AL Naghi company
AL Madinah road
opposite to Fetihi center,
Al Forsan
P.O. Box: 269
21411 Jeddah
☎ 02- 651 8670

Schweiz/Suisse/Svizzera

Telion AG, Rütistrasse 26,
8952 Schlieren,
☎ 0844-88 40 10

Singapore (Republic of)

Beste (S) Pte. Ltd.,
No. 6 Tagore Drive,
03-04 Tagore Industrial Building,
Singapore 787623,
☎ (65) 6552 2422

Slovakia

Techno Servis Bratislava,
Bajzova 11/A,
82108 Bratislava,
☎ (02) 55 56 37 49

Slovenia

Iskra Prins d.d.
Rozna dolina c. IX/6
1000 Ljubljana,
☎ 386 01 476 98 00

South Africa (Republic of)

Fixnet After Sales Service,
159 Queen Street, Kensington South,
P.O. Box 751770,
Johannesburg 2094,
Kensington South,
☎ 11 615 6765

St. Maarten

Rupchand Sons n.v. (ram's),
P.O. Box 79,
Philipsburg,
Netherlands Antilles,
☎ 52 29 31

St. Thomas

Boolchand's Ltd.,
31 Main Street,
P.O. Box 5667,
00803 St. Thomas,
US Virgin Islands,
☎ 340 776 0302

Suomi

Gillette Group Finland Oy,
P.O. Box 9,
Niittykatu 8, PL 9, 02200 Espoo,
☎ 09-45 28 71

Sverige

Gillette Group Sverige AB,
Dept. Sweden,
Stockholm Gillette
Räsundavägen 12,
Box 702,
16927 Solna,
☎ 020-21 33 21

Syria

Ahmed Hadaya Company
Hadaya building
Ain Keresh
Unisyria, P.O. Box 35002,
Damascus,
☎ 963 011-231433

Taiwan

Audio & Electr. Supplies Ltd.,
Brothers Bldg., 10th Floor,
85 Chung Shan N Rd., Sec. 1,
Taipei (104),
☎ (886) 02 2523 3283

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Yemen (Republic of)

Saba Stores for Trading,
26th September Street,
P.O. Box 5278,
Taiz,
☎ 4-25 23 88

Yugoslavia

BG Elektronik,
Bulevar kralja Aleksandra 34,
11000 Beograd,
☎ 11 3240 030

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Driva S.A.,
Gillette de Venezuela S.A.,