
L I F E I N B A L A N C E


# SOEHNLE

Downloaded from [www.vandenborre.be](http://www.vandenborre.be)

BEDIENUNGSANLEITUNG	2
USER MANUAL	4
MODE D'EMPLOI	6
ISTRUZIONI PER L'USO	8
GEBRUIKSAANWIJZING	10
INSTRUCCIONES DE MANEJO	12
MANUAL DE INSTRUÇÕES	14
BRUGSANVISNING	16
BRUKSANVISNING	18
INSTRUKCJA OBSŁUGI	20
NÁVOD NA OBSLUHU	22
NÁVOD K POUŽITÍ	24
ИНСТРУКЦИЯ ПО ЭКСПЛУАТАЦИИ	26
KULLANMA KILAVUZU	28
NAUDOJIMO INSTRUKCIJA	30
УПЪТВАНЕ ЗА УПОТРЕБА	32
دليل الاستخدام	34

Shape Sense Profi 300

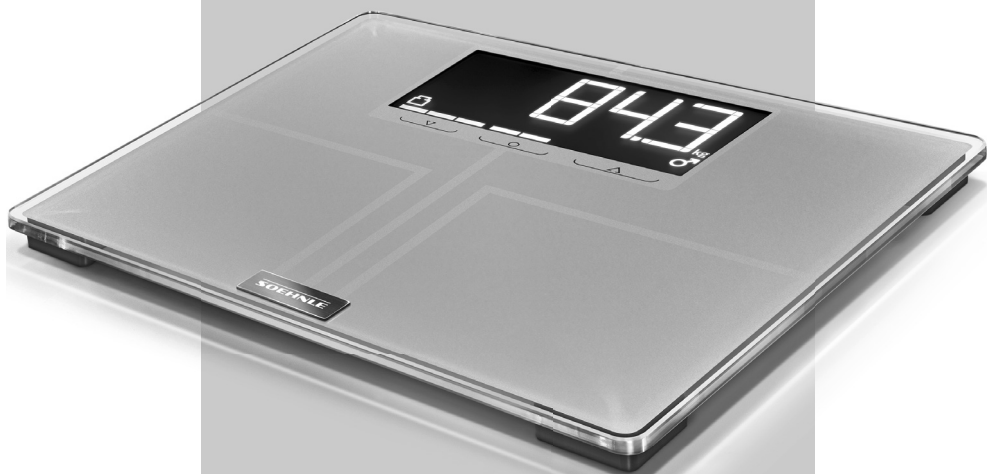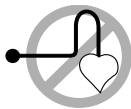- Nicht geeignet für Personen mit elektronischen Implantaten (z.B. Herzschrittmacher).
- It is not suitable for persons with electronic implants (pacemaker, etc.).
- Les personnes qui portent des implants médicaux (p. ex. un stimulateur cardiaque etc.) ne doivent pas utiliser cet appareil.
- Non è adatto per persone con impianti elettronici (pacemaker, ecc.).

- Niet geschikt voor personen met elektronische implantaten (pacemaker, enz.).
- No es apropiado para personas con implantes electrónicos (marcapasos, etc.).
- Não é adequado para pessoas com implantes eletrônicos (pacemakers, etc.).
- Ej lämplig för personer med elektroniska implantat (pace-maker etc).
- Ikke egnet til folk med elektroniske implantater (pacemakere ect).
- Nie nadaje się dla osób, którym wszczepiono elektroniczne implanty (np. rozruszniki serca).
- Prístroj nie je vhodný pre osoby s elektronickými implantátmi (napr. kardiostimulátorom).
- Není vhodný pro osoby s elektronickými implantáty (kardiostimulátory, atd.).
- Никогда не используйте данный прибор вместе с электронными медицинскими приборами, (например кардиостимулятор, и. т. п.)
- Elektronik parçalar (kalp pili, vs.) taşıyan insanlar için uygun değildir.
- Netinka elektroninius implantus (pvz. širdies stimulatorius) turintiems žmonėms.
- Уредът не е подходящ за хора с електронни импланти (напр. сърдечни пейсмейкъри).
- غير مناسب للأشخاص الذين يستخدمون أجهزة إلكترونية مزروعة (مُنظّم ضربات القلب).

Abweichung vom Normalwert • Deviation from the normal value • Déviation de la valeur normale • Deviazione dal valore normale • Afwijking van de normale waarde • La desviación del valor normal • Desvio do valor normal • Afvigelse fra den normale værdi • Avvikelse från normalvärdet • Odchyłka od běžné hodnoty • Odchyłka od normalnej hodnoty • Odchylenie od wartości normalnej • Отклонение от нормального значения • Normal degerden sapma • Nuokrypsiu nuo normaliosios vertės • Отклонение от нормальной стоимости • انحراف عن القيمة العادية

	1	2	3	4	5	6	7	8	9
	-			☺			+		
<b>BMI</b>	< 17,5	17,5	18	18,5	20,5	23	25,1	27,5	> 30
☹☹	< -12 %	- 12 %	- 7 %	-3 %	0 %	+ 3 %	+ 7 %	+ 12 %	> + 12 %
💧	< - 6 %	- 6 %	- 4 %	-2 %	0 %	+ 2 %	+ 4 %	+ 6 %	> + 6 %
🧑	< - 6 %	- 6 %	- 4 %	-2 %	0 %	+ 2 %	+ 4 %	+ 6 %	> + 6 %

Alter/Age	♀		☹☹		💧🧑		♂		☹☹		💧🧑	
	-	☺	+	++	☺	☺	-	☺	+	++	☺	☺
16-18	< 15	15-25	25-33	> 33	> 52,5	> 28	< 8	8-18	18-24	> 24	> 58,5	> 38
18-30	< 20	20-29	29-36	> 36	> 49,5	> 28	< 8	8-18	18-24	> 24	> 57	> 38
30-40	< 22	22-31	31-38	> 38	> 47,5	> 28	< 11	11-20	20-26	> 26	> 56	> 38
40-50	< 24	24-33	33-40	> 40	> 46,5	> 28	< 13	13-22	22-28	> 28	> 55	> 38
50-60	< 26	26-35	35-42	> 42	> 46	> 28	< 15	15-24	24-30	> 30	> 54	> 38
60+	< 28	28-37	37-47	> 47	> 45	> 28	< 17	17-26	26-34	> 34	> 53	> 38

Körperfettanteil, body fat content, taux de graisse, percentuale di grasso, lichaamsvetaandeel, porcentaje de la grasa corporal, percentagem de gordura corporal, fedtprocent, andelen kroppsfett, test zsrirtartalmát, zawartość tłuszczu w ciele, podíl tělesného tuku, доля внутреннего жира, үцүттәки yağ oranını, Kūno riebalų procentinė dalis, Процент на телесната машина, نسبة الدهون في الجسم


Körperwasseranteil, body water content, taux d'eau dans le corps, percentuale acqua corporei, lichaamswateraandeel, porcentaje del agua corporal, percentagem de água corporal, vandprocent, andelen kroppsvatten, test vızartalmát, zawartość wody w ciele, podíl vody v těle, доля тканевой жидкости, үцүттәки су oranını, K no vandens procentinis, Процент на телесната вода, نسبة الماء في الجسم

Muskelanteil, muscle proportion, pourcentage de muscles, massa muscular, spieraandeel, porción de los músculos, percentagem de músculo, muskler andelen, muskellandel, izom arány, masy mięśniowej, podíl svalů, доля мышечной, үцүттәки kas oranını, raumen procentin dalis, Процент на мускулите, نسبة العضلات

# A
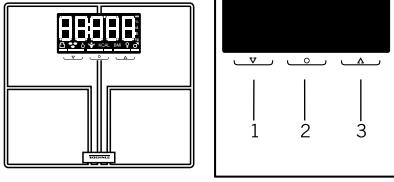
# B

- 
- 
- 
- 

# C

- 
- 
- 
- 
- 
- 
- 

# D

- 
- 
- 
- 
- 
- 
- 

# E

- - 
  -
- New batteries  
3 x 1,5 V AAA

## A. Bedienelemente

1. Minus (▼)
2. Bestätigen (○)
3. Plus (▲)

**Sensor-Touch:**

**Tasten nur sanft berühren!**

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1. Batterien (3 x 1,5 V AAA) einlegen.

### Hinweis:

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1. Insert batteries.

### Note:

Place the body analysis scale immediately on a solid and even surface (as long as „0,0“ is still displayed) and wait until the scale turns off automatically.

Only then begin with data input. Note, it might occur that the initial scale weight is being added erroneously to the personal weight process.

The same may happen during regular operation, if scale is held in the hand prior to the weight process.

2. Place scale level on a hard surface (not carpet). This applies for all measurements.
3. Cleaning and Care: use only a lightly damp cloth for cleaning. Do not use solvents or abrasives. Never submerge scale in water.
4. A country-specific option, to switch from kg/cm to lb/in or st/in, is available on the back of the scale by pushing key.

**Attention!** Slipping hazard in case of wet surface.

## C. Data input

1. Push ○ key.
2. Use ▲/▼ keys to select memory slot. Use ○ to confirm. → P
3. Then use ▲/▼ keys enter body height and confirm. → H
4. Enter the age and confirm. → A
5. Select gender and confirm. → S
6. To finish, set activity level and confirm (Press ○ key 3 sec.). → AC

1	Hardly any physical activity	
2	Slight, active, sitting and standing activity	
3	House and yard work, mainly standing, sitting only occasionally.	
4	Athletes and people with plenty of movement	intensive training: 5 - 10 h/week
5	Heavy manual labourer, highly active athletes	intensive training: min. 10 h/week

**Athlete Mode:** The respective calorie requirements can be precisely determined by setting personal activity level.

7. Place the body analysis scale immediately on the floor, wait until screen displays „0,0“. Step barefoot on the scale.

The initial measurement is indispensable and serves to store data for later automatic user recognition. If user does not step on the scale, data cannot be saved and must be entered again.

The first body analysis is performed after a user steps on the scale.

This scale is equipped with a Timeout mode. It will automatically switch off after approx. 40 seconds, if no keys are pressed.

Data must be entered again after battery change.

## D. Body analysis

1. Place scale on solid and even surface and wait until scale switches off automatically. Step barefoot on the scale. Stand still.
2. The determined body weight is displayed. The automatic recognition of user follows.
3. If the determined values cannot be clearly allocated to one user, the possible memory slots will be displayed.

If the determined weight can be allocated to more than 2 users, other possible memory slots will be displayed after 3 seconds.

Use the ▲/▼ keys (▼= left, ▲ = right) to confirm the memory slot displayed there.

The scale measures body values (body fat, body water, muscle mass), basic metabolic rate (kcal) and BMI (Body Mass Index).

The 9-segment status bar shows the respective deviation from the normal value (see overleaf).

The scale turns off automatically.

## Classification of body weight according to BMI

(Source: World Health Organisation, 2016)

BMI (kg/m <sup>2</sup> )	classification
< 18.5	Underweight
18.5 - 24.9	Normal weight
25.0 - 29.9	Pre-obesity
30.0 - 34.9	Obesity class I
35.0 - 39.9	Obesity class II
> 40.0	Obesity class III

## E. Notifications

1. Replace batteries
2. Overload: max. 200 kg
3. The measured value is higher (H) or lower (L) than the limit.

	L  H	L  H	L  H
♀	< 10 % > 60 %	< 5 % > 90 %	< 5 % > 90 %
♂	< 4 % > 60 %	< 5 % > 90 %	< 5 % > 90 %

## F. Technical data

Max. 200 kg / 100 g

Max. 440 lb / 0.2 lb

Max. 31 st / 0.2 lb

Battery requirement: 3 x 1,5 V AAA

## Guarantee

Leifheit AG grants you 3 years guarantee as of the date of purchase (or for order upon receipt of the goods) for quality product at hand. The guarantee claims must be asserted immediately upon occurrence of defect within the guarantee period. The guarantee is limited to quality of the products.

The guarantee does not cover:

- (1) wear caused by proper use or natural causes,
- (2) damage caused by improper use or handling (e.g. impact, bump, fall),
- (3) damage caused by non-compliance with stated operating instructions,
- (4) the rechargeable battery.

In case of a guarantee claim, Leifheit will at its own discretion, provide either repair of defective parts or replace the product. If a repair is not possible and if an identical product is no longer available within the product range for replacement, you will receive a replacement product that is,

insofar as possible, of equal value. In case of a guarantee claim, the purchase price cannot be reimbursed. This guarantee also does not grant compensation rights. Furthermore, this guarantee does not grant any claims for damages.

To assert the guarantee, please contact seller where you purchased the product and present the defective product, purchase receipt (copy) and guarantee slip.

Guarantee is valid globally.

Your statutory rights, in particular warranty rights, continue to be valid and are not limited by this warranty.

## EC Conformity

This device complies with the applicable EC Directive 2014/30/EC.

You will find the EU compliance declaration on [www.soehnle.com](http://www.soehnle.com).

## Battery disposal EC Directive 2008/12/EC

Batteries are not a part of your regular household waste. You must return batteries to your municipality's public collection points or wherever batteries of respective type are being sold.

Batteries contain hazardous materials are marked with these symbols:

Pb = battery contains lead

Cd = battery contains cadmium

Hg = battery contains mercury

## Disposal of electric and electronic devices EC Directive 2002/96/EC

This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

## Consumer service

If you have questions or comments, the following points of contact are available:

<b>Germany</b>	Ph.: (08 00) 5 34 34 34
<b>International</b>	Ph.: +49 26 04 97 70
<b>Monday to Friday</b>	08:30 am to 12:00 am

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LEIFHEIT AG