L I F E I N B A L A N C E

SOEHNLE



BEDIENUNGSANLEITUNG	2
USER MANUAL	4
MODE D'EMPLOI	6
ISTRUZIONI PER L'USO	8
GEBRUIKSAANWIJZING	10
INSTRUCCIONES DE MANEJO	12
MANUAL DE INSTRUÇÕES	14
BRUGSANVISNING	16
BRUKSANVISNING	18
INSTRUKCJA OBSŁUGI	20
NÁVOD NA OBSLUHU	22
NÁVOD K POUŽITÍ	24
ИНСТРУКЦИЯ ПО ЭКСПЛУАТАЦИИ	26
KULLANMA KILAVUZU	28
СОБСТВЕНИЦИТЕ НА ПОТРЕБИТЕЛЯ	30
SAVININKAI VADOVAS	32
دليل مالكي	34

Style Sense Multi 200







Deutsch	2	Poertuguês	14	Русский	26
English	4	Dansk	16	Türkçe	28
Français	6	Svenska	18	Lietuvos	30
Italiano	8	Český	20	български	32
Nederlands	10	Slowensky	22	اللغة العربية	34
Español	12	Polsku	24		






A


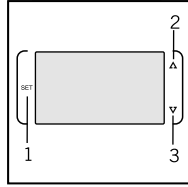
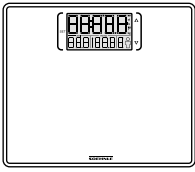
B

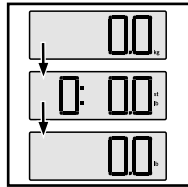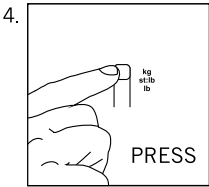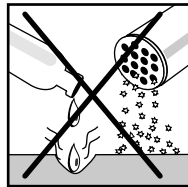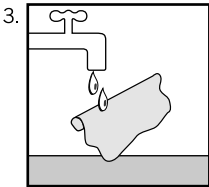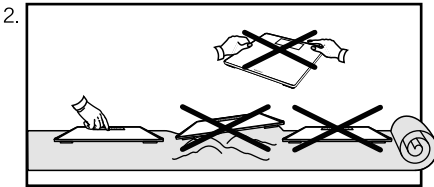
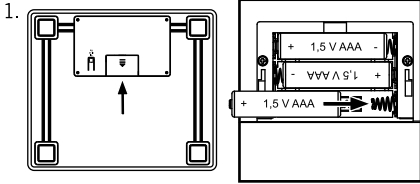
C

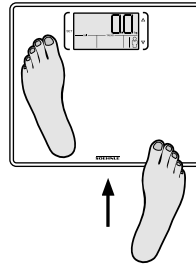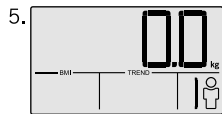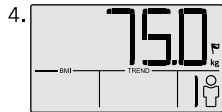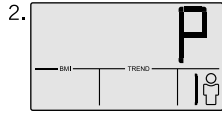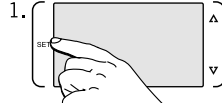D

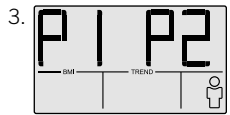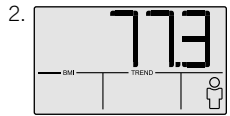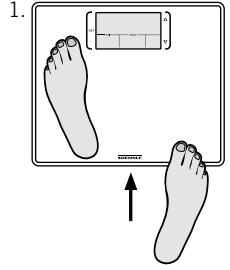E

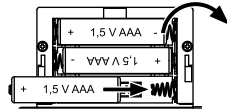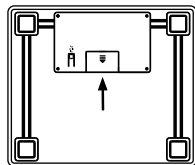New batteries
3 x 1,5 V Size AAA

A. Bedienelemente

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1. Insert batteries.

Note:

Place the body analysis scale immediately on solid and even surface (as long as „0,0“ is still displayed) and wait until the scale turns off automatically.

Only then begin with data input. Note, it might occur that the initial scale weight is being added erroneously to the personal weight process.

The same may happen during regular operation, if scale is held in the hand prior to the weight process.

2. Place scale on solid and even surface (not carpet). This applies for all measurements.
3. Cleaning and Care: use only a lightly damp cloth for cleaning. Do not use solvents or abrasives. Never submerge scale in water.
4. A country-specific option, to switch from kg/cm to lb/in or st/in, is available on the back of the scale by pushing key.

Attention! Slipping hazard in case of wet surface.

C. Data input

1. Push SET key.
2. Use Δ/∇ keys to select memory slot and confirm (SET).
3. Then use Δ/∇ keys to enter body height and confirm.
4. Set target weight and confirm (Press SET key 3 sec.).
5. Place the body analysis scale immediately on the floor, wait until screen displays „0,0“ and then step on the scale.

The initial measurement is indispensable and serves to store data for later automatic recognition. If no user steps on the scale, data input must be repeated.

The first body analysis is performed after a user steps on the scale.

This scale is equipped with a timeout mode. It will automatically turn off after approx. 40 seconds, if no keys are pressed.

Data must be entered again after battery change.

D. Body analysis

1. Place scale on solid and even surface and wait until scale switches off automatically. Step on scale. Please stand still.
2. Determined body weight is displayed. Automatic user recognition follows.
3. If the determined values cannot be clearly allocated to one user, the possible memory slots will be displayed.
If the determined weight can be allocated to more than 2 users, other possible memory slots will be displayed after 3 seconds.
Use the Δ/∇ keys (Δ = left, ∇ = right) to confirm the memory slot displayed there.
4. The determined body weight, BMI (Body Mass Index) and the difference to the target weight is displayed.
The scale turns off automatically.

Classification of body weight according to BMI

(Source: World Health Organisation, 2016)

BMI (kg/m ²)	classification
< 18.5	Underweight
18.5 - 24.9	Normal weight
25.0 - 29.9	Pre-obesity
30.0 - 34.9	Obesity level I
35.0 - 39.9	Obesity level II
> 40.0	Obesity level III

E. Notifications

1. Replace batteries
2. Overload: max. 180 kg

F. Technical data

Max. 180 kg / 100 g

Max. 396 lb / 0.2 lb

Max. 28 st / 0.2 lb

Battery requirement: 3 x 1,5 V AAA

Guarantee

Leifheit AG grants you 5 years guarantee as of the date of purchase (or for order upon receipt of the goods) for the quality product at hand. The guarantee claims must be asserted immediately upon occurrence of the defect within the guarantee period.

The guarantee is limited to the quality of the products.

The following are excluded from the guarantee:

- (1) wear caused by proper use or natural causes,
- (2) damage caused by improper use or handling (e.g. impact, bump, fall),
- (3) damage caused by non-compliance with stated operating instructions,
- (4) rechargeable battery.

In case of a guarantee claim, Leifheit will at its own discretion provide either repair of defective parts or replace the product. If a repair is not possible and if an identical product is no longer available within the product range for replacement, you will receive a replacement product that is, insofar as possible, of equal value. In case of a guarantee claim, the purchase price cannot be reimbursed. This guarantee also does not grant compensation rights. Furthermore, this guarantee does not grant any claims for damages.

To assert the guarantee, please contact the seller where you purchased the product and present defective product, purchase receipt (copy) and guarantee slip.

The guarantee is valid globally.

Your statutory rights, in particular warranty rights, continue to be valid and are not limited by this warranty.

EC Conformity

 This device complies with the applicable EC Directive 2014/30/EC.

You will find the EU compliance declaration on www.soehnle.com.

Battery disposal

EC Directive 2008/12/EC

Batteries are not a part of your regular household waste. You must return batteries to your municipality's public collection points or wherever batteries of respective type are being sold.

Batteries containing hazardous materials are marked with these symbols:

Pb = battery contains lead

Cd = battery contains cadmium

Hg = battery contains mercury

Disposal of electric and electronic devices EC Directive 2002/96/EC

This product is not to be treated as regular household waste but must be returned to  a collection point for recycling electric and electronic devices. Further information is provided by municipality, municipality's waste disposal services, or retailer where product was purchased.

Consumer service

If you have questions or comments, the following points of contact are available:

Germany	Ph.: (08 00) 5 34 34 34
International	Ph.: +49 26 04 97 70
Monday to Friday	08:30 am to 12:00 am

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LEIFHEIT AG